

# Series 2000XL


## 3-sectioned multipurpose table with electric motorized height adjustment

... for examination, rectoscopy, gynaecology ...

Model 2645-00XLE white powder coated (RAL 9010)

Model 2645-03XLE white aluminium powder coated (RAL 9006)

Model 2645-04XLE grey aluminium powder coated (RAL 9007)

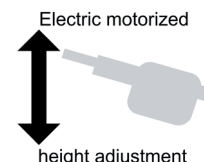

Model 2645-00XLE  
with option :  
Movable on twin castors  
with central lock (1180N);  
upper picture additionally:  
Leg supports (1510);  
Upholstery colour:  
64-173 smaragd



### Basic equipment:

- Movable on comfort twin castors, Ø 100 mm, each with lock
- Hand grips on both sides of the seat section
- Drawable basin (stainless steel) under seat, basin support: HPL panel
- Including clamps for leg supports (leg supports = available as an option)
- Electromotive drive control by foot switch (option: hand switch)
- Integrated safety function for automatic blocking of electromotive drive (art. no. 1999, mandatory equipment)
- Freewheeling downwards, automatic declutch in case of resistance
- Dewert-hygiene-equipment



### Upholstery:

Lying surface 3-teilig:

- Head part 30° positive adjustable by 2 metal rachets
- Length head part: 550 mm
- Middle part 50° positive adjustable by rack system
- Length middle part: 930 mm
- Seat section with half-round cut-out, manually inclinable (20°)
- Length seat section: 450 mm
- Thickness of upholstery: 80 mm

### Technical data:

- Electronical height adjustment: 620 mm to 1060 mm in 22 sec. by AC-motor with a pressure capacity of 10000 N
- Width of lying surface: 650 mm
- Length of lying surface: 1980 mm
- Load capacity: 220 kg; according to IEC 60601 tested with fourfold static load



### Design:

2 different variants for artificial leather eligible for upholstery covering:

- Skai® Pandoria Plus with velvet-smooth surface
- Skai® Toronto EN with innovative staynu-technology



Available as an option (see pricelist for additional options)

### Mobility



Art. no. 1182  
Without basic mobility



Art. no. 1180N  
Movable on twin castors, Ø 100 mm with central lock, basic height +20 mm

### Additional motor control options



Art. no. 1920E  
Hand switch (low voltage), with magnet on the back to fasten



Art. no. 1910E  
Foot switch (low voltage)



Art. no. 1950E  
Foot switch rail, lateral, on both sides, security adjustment

### Side guards



Art. no. 1040  
1 side guard tiltable, chromium-plated, completely made of steel tube, length: 700 mm



Art. no. 1005  
1 side guard to fit on side rail, with attachment clamp, chromium-plated, length: 800 mm

### Add-on elements



Art. no. 1410  
1 piece of side rail, length: 350 mm



Art. no. 1430  
1 attachment clamp for side rail, suitable for pipe Ø 18 mm

### Leg supports



Art. no. 1510  
1 pair of leg supports



Art. no. 1520  
1 pair of leg supports according to Goepel

### Nose opening



Art. no. 1610  
Nose opening, in head part

Art. no. 1620  
Upholstery part to close



Art. no. 1940E  
1 foot switch fixation mounted at underframe

### Paper roll holder



Art. no. 1265K  
Paper roll holder, mounted under head part (not protruding), for width of lying surface 700 mm

### Head part options



Art. no. 1655  
Head part stepless adjustable by gas spring from 0° to +60°



Art. no. 1450  
1 Infusion pole for side rail, without attachment clamp

### Upholstery options

Art. no. 1770  
Soft covering for pleasant resting comfort, for width 700 mm, thickness 80 mm

Art. no. 1772  
Upholstery medium-soft (cold foam cover), for width 700 mm, thickness 70 mm

