

Tables for echocardiography

2-sectioned table for echocardiography with electric motorized height adjustment ... for examination, treatment ...

Model 2550-00XLE and 2551-00XLE white powder coated (RAL 9010)
 Model 2550-03XLE and 2551-03XLE white aluminium powder coated (RAL 9006)
 Model 2550-04XLE and 2551-04XLE grey aluminium powder coated (RAL 9007)



2 upper pictures:
 Model 2550-00XLE with options: Movable on twin castors with central lock (1180N), paper roll holder (1270K); upholstery colour 64-717 deep sea



Basic equipment:

- Table not movable, including height adjustment screws to compensate for uneven ground
- In head part: Cut-out for echocardiography and upholstery part to close (in patient supine position left)
- Electromotive drive control by foot switch (option: hand switch)
- Integrated safety function for automatic blocking of electromotive drive
- Freewheeling downwards, automatic declutch in case of resistance
- Plug connector for potential equalization in acc. to DIN 42801
- Dewert-hygiene-equipment



Model 2550-00XLE with options: Movable on comfort twin castors, each with lock (1120N-10), storage basket (1850), electric motorized head part adjustment (1665-10)

Upholstery:

- Lying surface 2-sectioned:
- Head part stepless upwards (+50°) adjustable by gas spring
 - Length head part: 800 mm
 - In head part: Cut-out for echocardiography and upholstery part to close, in patient supine position left, measures: L. 300 mm x W. 250 mm
 - Fixed foot part
 - Length foot part: 1150 mm
 - Upholstery thickness: 60 mm


Technical data:

- Electronical height adjustment: 500 mm to 940 mm in 22 sec. by AC-motor with a pressure capacity of 10000 N
- Width of lying surface: Model 2550-..XLE: 700 mm | Model 2551-..XLE: 800 mm
- Length of lying surface: 1950 mm
- Load capacity: 250 kg; according to IEC 60601 tested with fourfold static load



Design:

- 2 different variants for artificial leather eligible for upholstery covering:
- Skai® Pandoria Plus with velvet-smooth surface
 - Skai® Toronto EN with innovative staynu-technology

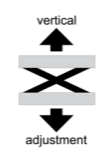

Available as an option (see pricelist for additional options)

Mobility



Art. no. 1120N-10
 Movable on comfort twin castors Ø 100 mm, each with total lock, basic height -15 mm



Art. no. 1180N
 Movable on twin castors, Ø 100 mm with central lock, basic height +20 mm



Art. no. 1185
 Movable on comfort twin castors, Ø 125 mm with central lock, basic height + 50 mm



Art. no. 1045N
 Wheel-lifting/fixing mechanism with 4 foot levers for controlled lowering

Additional motor control options



Art. no. 1920E
 Hand switch (low voltage), with magnet on the back to fasten



Art. no. 1910E
 Foot switch (low voltage)



Art. no. 1940E
 1 foot switch fixation mounted at underframe



Art. no. 1950E
 Foot switch rail, lateral, on both sides, security adjustment




Art. no. 1170
 1 wall protection bumper Ø 150 mm

Head part options



Art. no. 1665-10
 Electric motorized head part adjustment to +60°

Paper roll holder

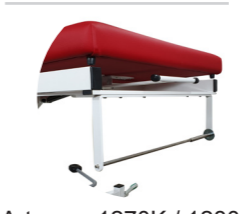

Art. no. 1270K / 1280K
 Paper roll holder, mounted under head part (not protruding), for width 700 mm / 800 mm



Art. no. 1270F / 1280F
 Paper roll holder, mounted at foot part (protruding: 70 mm), for width 700 mm / 800 mm

Upholstery options

Art. no. 1770 / 1780
 Soft covering for pleasant resting comfort, for width 700 mm / 800 mm, thickness 80 mm

Art. no. 1772 / 1782
 Upholstery medium-soft (cold foam cover), for width 700 mm / 800 mm, thickness 70 mm

Add-on elements



Art. no. 1410
 1 piece of side rail, length: 350 mm



Art. no. 1850
 Storage basket

Side guards



Art. no. 1080-50
 1 side guard lateral lowerable, chromium-plated, completely made of steel tube, length: 1060 mm, length on cut-out side: 820 mm



Art. no. 1040
 1 side guard tiltable, chromium-plated, completely made of steel tube, length: 700 mm



Art. no. 1005
 1 side guard to fit on side rail, with attachment clamp, chromium-plated, length: 800 mm

