


Series 2800XL

3-sectioned therapy table with electric motorized height adjustment

Model 2815-00XLE white powder coated (RAL 9010)

Model 2815-03XLE white aluminium powder coated (RAL 9006)

Model 2815-04XLE grey aluminium powder coated (RAL 9007)

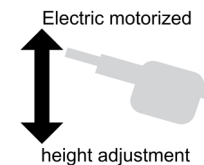

without
picture



Model 2815-00XL/H (hydraulic height adjustment) with option:
Movable on twin castors with central lock (046);
upholstery colour: 64-716 black

Basic equipment:

- Lying surface 3-sectioned
- Foot part stepless upwards (+75°) adjustable to sitting position
- Table not movable, including height adjustment screws to compensate for uneven ground
- Electromotive drive control by foot switch (option: hand switch)
- Integrated safety function for automatic blocking of electromotive drive
- Freewheeling downwards, automatic declutch in case of resistance
- Dewert-hygiene-equipment



Under part in case
of electric motorized
height adjustment

Upholstery:

Lying surface 3-sectioned:

- Head part stepless up- (+45°) and downwards (-35°) adjustable by gas spring
- Length head part: 540 mm, with nose opening and upholstery part to close
- Fixed middle part
- Length middle part: 530 mm
- Foot part stepless upwards (+75°) adjustable to sitting position by gas spring
- Length foot part: 880 mm
- Thickness of upholstery: 60 mm

Technical data:

- Electronical height adjustment: 510 mm to 950 mm in 22 sec. by AC-motor with a pressure capacity of 10000 N
- Width of lying surface: 700 mm
- Length of lying surface: 1950 mm
- Load capacity: 250 kg; according to IEC 60601 tested with fourfold static load



Design:

2 different variants for artificial leather eligible for upholstery covering:

- Skai® Pandoria Plus with velvet-smooth surface
- Skai® Toronto EN with innovative staynu-technology



Available as an option (see pricelist for additional options)

Mobility



Art. no. 040N
Movable on comfort twin
castors Ø 100 mm,
each with total lock,
basic height -15 mm



Art. no. 046
Movable on twin castors,
Ø 100 mm with central
lock, basic height +20 mm



Art. no. B1185
Movable on comfort twin
castors, Ø 125 mm with
central lock,
basic height + 50 mm



Art. no. 045N
Wheel-lifting/fixing mechanism
with 4 foot levers for controlled
lowering

Additional motor control options



Art. no. 090E
Hand switch (low voltage),
with magnet on the back to fasten



Art. no. 091E
Foot switch (low voltage)



Art. no. 092E
1 foot switch fixation mounted at
underframe



Art. no. 094E
Foot switch rail, lateral, on both sides,
security adjustment:
Press the switch rail down = the table is lifted
Lift the switch rail up = the table is lowered



Art. no. 095E
Foot switch rail, surround,
security adjustment:
Press the switch rail down = the table is lifted
Lift the switch rail up = the table is lowered



Art. no. B1997E
Foot switch cover,
made of brushed VA-steel

Nose opening



Art. no. 1615-EF
Loss of nose opening and upholstery
part to close

Paper roll holder



Art. no. 010K
Paper roll holder, mounted under head
part (not protruding), for width 700 mm,
not in case of 3-sectioned head part
(Art. no. 056)



Art. no. 010K/F
Paper roll holder, mounted at foot part
(not protruding), for width 700 mm

Add-on elements



Art. no. B1410
1 piece of side rail, length:
350 mm



Art. no. B1430
1 attachment clamp for
side rail,
suitable for pipe: Ø 18 mm

Arm rests

without
picture

Art. no. 088
1 arm rest, fixed, swivelling
and removable, additional
revolving and tiltable



Art. no. B1440
1 arm rest for side rail, without
attachment clamp, swivelling and
removable, additional revolving and
tiltable

