


Series 2800XL

3-sectioned therapy table with hydraulic height adjustment

Model 2815-00XL/H white powder coated (RAL 9010)

Model 2815-03XL/H white aluminium powder coated (RAL 9006)

Model 2815-04XL/H grey aluminium powder coated (RAL 9007)

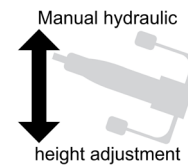

All pictures:
Model 2815-00XL/H
with options:
Movable on twin castors with
central lock (046), paper roll holder
(010K);
upholstery colour: 64-711 ashgrey

Basic equipment:

- Lying surface 3-sectioned
- Foot part stepless upwards (+75°) adjustable to sitting position
- Table not movable, including height adjustment screws to compensate for uneven ground
- Pedals for foot-controlled height positioning on both sides of the table
- Dewert-hygiene-equipment



Upholstery:

Lying surface 3-sectioned:

- Head part stepless up- (+45°) and downwards (-35°) adjustable by gas spring
- Length head part: 540 mm, with nose opening and upholstery part to close
- Fixed middle part
- Length middle part: 530 mm
- Foot part stepless upwards (+75°) adjustable to sitting position by gas spring
- Length foot part: 880 mm
- Thickness of upholstery: 60 mm



Technical data:

- Manual foot operated height adjustment: 530 mm to 950 mm by a hydraulic pump with a pressure capacity of 10000 N
- Width of lying surface: 700 mm
- Length of lying surface: 1950 mm
- Load capacity: 250 kg; according to IEC 60601 tested with fourfold static load



Design:

2 different variants for artificial leather eligible for upholstery covering:

- Skai® Pandoria Plus with velvet-smooth surface
- Skai® Toronto EN with innovative staynu-technology



Available as an option (see pricelist for additional options)

Mobility



Art. no. 040N
Movable on comfort twin
castors Ø 100 mm,
each with total lock,
basic height -15 mm



Art. no. 046
Movable on twin castors,
Ø 100 mm with central
lock, basic height +20 mm



Art. no. B1185
Movable on comfort twin
castors, Ø 125 mm with
central lock,
basic height + 50 mm



Art. no. 045N
Wheel-lifting/fixing mechanism
with 4 foot levers for controlled
lowering



Paper roll holder



Art. no. 010K
Paper roll holder, mounted under head
part (not protruding), for width 700 mm,
not in case of 3-sectioned head part
(Art. no. 056)



Art. no. 010K/F
Paper roll holder, mounted at foot part
(not protruding), for width 700 mm

Upholstery options

Art. no. 007
Width of upholstery 800 mm

Art. no. B1770
Soft covering for pleasant resting comfort,
for width 700 mm,
thickness 80 mm


Art. no. B1772
Upholstery medium-soft (cold foam
cover), for width 700 mm,
thickness 70 mm

Head part options



Art. no. 056
Head part 3-sectioned, arm rests
downwards adjustable from -25° to
+45° and removable
(art. no. 010K not possible)
measures:
Head part: L. 540 mm x W. 370 mm
Arm rests: L. 500 mm x W. 160 mm

Nose opening

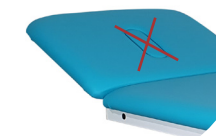

Art. no. 1615-EF
Loss of nose opening and upholstery
part to close

Add-on elements



Art. no. B1410
1 piece of side rail, length:
350 mm



Art. no. B1430
1 attachment clamp for
side rail, suitable for pipe:
Ø 18 mm

Arm rests

without
picture

Art. no. 088
1 arm rest, fixed,
swivelling and
removable, additional
revolving and tiltable



Art. no. B1440
1 arm rest for side rail, without
attachment clamp, swivelling
and removable, additional
revolving and tiltable