


Series 2800XL

3-sectioned therapy table with electric motorized height adjustment

Model 2820-00XLE white powder coated (RAL 9010)

Model 2820-03XLE white aluminium powder coated (RAL 9006)

Model 2820-04XLE grey aluminium powder coated (RAL 9007)

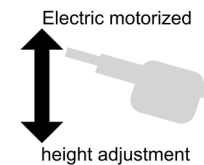

All pictures:
Model 2820-04XLE
with options:
Movable on twin castors with
central lock (046), side rail
(B1410);
upholstery colour: 64-716 black



Basic equipment:

- Lying surface 3-sectioned
- Head part can be hinged vertical
- Foot part stepless upwards adjustable to sitting position (+80°)
- Table not movable, including height adjustment screws to compensate for uneven ground
- Electromotive drive control by foot switch (option: hand switch)
- Integrated safety function for automatic blocking of electromotive drive
- Freewheeling downwards, automatic declutch in case of resistance
- Dewert-hygiene-equipment



Upholstery:

Lying surface 3-sectioned:

- Head part stepless up- (+30°) and downwards (-90°) adjustable by gas spring
- Length head part: 450 mm, with nose opening and upholstery part to close
- Fixed middle part
- Length middle part: 530 mm
- Foot part stepless upwards (+80°) adjustable to sitting position by gas spring
- Length foot part: 1020 mm
- Thickness of upholstery: 60 mm



Technical data:


- Electronical height adjustment: 510 mm to 950 mm in 22 sec. by AC-motor with a pressure capacity of 10000 N
- Width of lying surface: 700 mm
- Length of lying surface: 2000 mm
- Load capacity: 250 kg; according to IEC 60601 tested with fourfold static load



Design:

2 different variants for artificial leather eligible for upholstery covering:

- Skai® Pandoria Plus with velvet-smooth surface
- Skai® Toronto EN with innovative staynu-technology

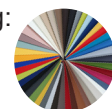

Available as an option (see pricelist for additional options)

Mobility



Art. no. 040N
Movable on comfort twin
castors Ø 100 mm,
each with total lock,
basic height -15 mm



Art. no. 046
Movable on twin castors,
Ø 100 mm with central
lock, basic height +20 mm



Art. no. B1185
Movable on comfort twin
castors, Ø 125 mm with
central lock,
basic height + 50 mm



Art. no. 045N
Wheel-lifting/fixing mechanism
with 4 foot levers for controlled
lowering



Additional motor control options



Art. no. 090E
Hand switch (low voltage),
with magnet on the back to fasten



Art. no. 091E
Foot switch (low voltage)



Art. no. 092E
1 foot switch fixation mounted at
underframe

Additional safety options



Art. no. 094E
Foot switch rail, lateral, on both sides,
security adjustment:
Press the switch rail down = the table is lifted
Lift the switch rail up = the table is lowered



Art. no. B1997E
Foot switch cover,
made of brushed VA-steel



Art. no. 1615-EF
Loss of nose opening and
upholstery part to close

Paper roll holder



Art. no. 010K
Paper roll holder, mounted under head
part (not protruding), for width 700 mm,
not in case of 3-sectioned head part
(Art. no. 056)



Art. no. 010K/F
Paper roll holder, mounted at foot part
(not protruding), for width 700 mm

Add-on elements



Art. no. B1410
1 piece of side rail, length:
350 mm



Art. no. B1430
1 attachment clamp for side rail,
suitable for pipe: Ø 18 mm



Art. no. B1440
1 arm rest for side rail, without
attachment clamp, swivelling and
removable, additional revolving and
tiltable