


Series 2800XL

3-sectioned therapy table with hydraulic height adjustment

Model 2820-00XL/H white powder coated (RAL 9010)

Model 2820-03XL/H white aluminium powder coated (RAL 9006)

Model 2820-04XL/H grey aluminium powder coated (RAL 9007)

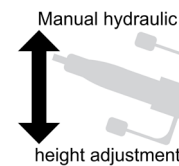

All pictures:
Model 2820-04XL/H
with options:
Movable on twin castors with
central lock (046), paper roll
holder (010K/F);
upholstery colour: 64-713 corn

Basic equipment:

- Lying surface 3-sectioned
- Head part can be hinged vertical
- Foot part stepless upwards adjustable to sitting position (+80°)
- Table not movable, including height adjustment screws to compensate for uneven ground
- Pedals for foot-controlled height positioning on both sides of the table
- Dewert-hygiene-equipment



Upholstery:

Lying surface 3-sectioned:

- Head part stepless up- (+30°) and downwards (-90°) adjustable by gas spring
- Length head part: 450 mm, with nose opening and upholstery part to close
- Fixed middle part
- Length middle part: 530 mm
- Foot part stepless upwards (+80°) adjustable to sitting position by gas spring
- Length foot part: 1020 mm
- Thickness of upholstery: 60 mm

Technical data:

- Manual foot operated height adjustment: 530 mm to 950 mm by a hydraulic pump, pressure capacity of 10000 N
- Width of lying surface: 700 mm
- Length of lying surface: 2000 mm
- Load capacity: 250 kg; according to IEC 60601 tested with fourfold static load



Design:

2 different variants for artificial leather eligible for upholstery covering:

- Skai® Pandoria Plus with velvet-smooth surface
- Skai® Toronto EN with innovative staynu-technology



Available as an option (see pricelist for additional options)

Mobility



Art. no. 040N
Movable on comfort twin
castors Ø 100 mm,
each with total lock,
basic height -15 mm



Art. no. 046
Movable on twin castors,
Ø 100 mm with central
lock, basic height +20 mm



Art. no. B1185
Movable on comfort twin
castors, Ø 125 mm with
central lock,
basic height + 50 mm



Art. no. 045N
Wheel-lifting/fixing mechanism
with 4 foot levers for controlled
lowering



Paper roll holder



Art. no. 010K
Paper roll holder, mounted under head
part (not protruding), for width 700 mm



Art. no. 010K/F
Paper roll holder, mounted at foot part
(not protruding), for width 700 mm

Add-on elements



Art. no. B1410
1 piece of side rail, length:
350 mm




Art. no. B1430
1 attachment clamp for
side rail,
suitable for pipe: Ø 18 mm

Arm rests



Art. no. B1440
1 arm rest for side rail,
without attachment
clamp, swivelling and
removable, additional
revolving and tiltable

Nose opening

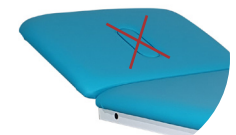

Art. no. 1615-EF
Loss of nose opening and upholstery
part to close