


Series 2800XL

4-sectioned therapy table with electric motorized height adjustment

Model 2895-00XLE white powder coated (RAL 9010)

Model 2895-03XLE white aluminium powder coated (RAL 9006)

Model 2895-04XLE grey aluminium powder coated (RAL 9007)

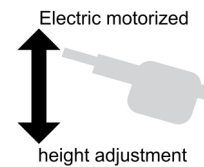

Model 2895-00XLE
with options:
Movable on twin castors with
central lock (046), paper roll
holder (010F/K);
upholstery colour: 64-153 platin

Basic equipment:

- Lying surface 4-sectioned
- Multifunctional and independent adjustment of upholstery sections
- Possibilities of positioning the surface parts: Roof shape, Trendelenburg, reverse Trendelenburg
- Table not movable, including height adjustment screws to compensate for uneven ground
- Electromotive drive control by foot switch (option: hand switch)
- Integrated safety function for automatic blocking of electromotive drive
- Freewheeling downwards, automatic declutch in case of resistance
- Dewert-hygiene-equipment



Upholstery:

Lying surface 4-sectioned:

- Head part stepless up- (+40°) and downwards (-35°) adjustable by gas spring
- Length head part: 420 mm; with nose opening and upholstery part to close
- Middle part stepless up- (+45°) and downwards (-15°) adjustable by gas spring
- Length middle part: 550 mm
- Foot part 2-sectioned, split transversal:
 - Upper foot segment stepless up- (+50°) and downwards (-15°) adjustable by gas spring
 - Lower foot segment stepless up- (+35°) and downwards (-30°) adjustable by gas spring
- Length foot part: Upper foot segment: 550 mm | Lower foot segment: 420 mm
- Thickness of upholstery: 60 mm



Technical data:


- Electronical height adjustment: 570 mm to 1010 mm in 22 sec.
by AC-motor with a pressure capacity of 10000 N
- Width of lying surface: 700 mm
- Length of lying surface: 1950 mm
- Load capacity: 200 kg; according to IEC 60601 tested with fourfold static load



Design:

2 different variants for artificial leather eligible for upholstery covering:

- Skai® Pandoria Plus with velvet-smooth surface
- Skai® Toronto EN with innovative staynu-technology

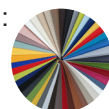


Available as an option (see pricelist for additional options)

Mobility



Art. no. 077N
Movable on twin castors
Ø 75 mm, each with total lock

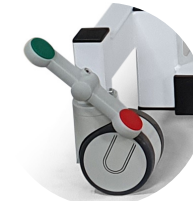

Art. no. 049
Movable on twin castors, Ø 100 mm
with central lock



Art. no. 044N
Wheel-lifting/fixing mechanism
with 4 foot levers for controlled
lowering

Additional motor control options



Art. no. 090E
Hand switch (low voltage),
with magnet on the back to fasten



Art. no. 091E
Foot switch (low voltage)



Art. no. 092E
1 foot switch fixation mounted at
underframe

Additional safety options



Art. no. B1997E
Foot switch cover,
made of brushed VA-steel

Paper roll holder



Art. no. 010F/K
Paper roll holder, mounted under head
part (not protruding), for width 700 mm

Nose opening



Art. no. 1615-EF
Loss of nose opening and upholstery
part to close

Upholstery options

Art. no. B1770
Soft covering for pleasant resting comfort,
for width 700 mm, thickness 80 mm

Art. no. B1772
Upholstery medium-soft (cold foam cover),
for width 700 mm, thickness 70 mm

