


Series 2800XL

4-sectioned therapy table with hydraulic height adjustment

Model 2895-00XL/H white powder coated (RAL 9010)

Model 2895-03XL/H white aluminium powder coated (RAL 9006)

Model 2895-04XL/H grey aluminium powder coated (RAL 9007)

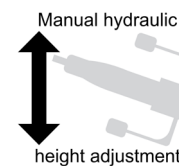

Under part in case of hydraulic height adjustment, with options:
Movable on twin castors with central lock (046)

Model 2895-00XLE (electric motorized height adjustment) with options:
Movable on twin castors with central lock (046),
paper roll holder (010F/K);
upholstery colour: 64-153 platin

without picture

Basic equipment:

- Lying surface 4-sectioned
- Multifunctional and independent adjustment of upholstery sections
- Possibilities of positioning the surface parts: Roof shape, Trendelenburg, reverse Trendelenburg
- Table not movable, including height adjustment screws to compensate for uneven ground
- Pedals for foot-controlled height positioning on both sides of the table
- Dewert-hygiene-equipment



Upholstery:

Lying surface 4-sectioned:

- Head part stepless up- (+40°) and downwards (-35°) adjustable by gas spring
- Length head part: 420 mm; with nose opening and upholstery part to close
- Middle part stepless up- (+45°) and downwards (-15°) adjustable by gas spring
- Length middle part: 550 mm
- Foot part 2-sectioned, split transversal:
 - Upper foot segment stepless up- (+50°) and downwards (-15°) adjustable by gas spring
 - Lower foot segment stepless up- (+35°) and downwards (-30°) adjustable by gas spring
- Length foot part: Upper foot segment: 550 mm | Lower foot segment: 420 mm
- Thickness of upholstery: 60 mm

Technical data:

- Manual foot operated height adjustment: 600 mm to 1020 mm by a hydraulic pump, pressure capacity of 10000 N
- Width of lying surface: 700 mm
- Length of lying surface: 1950 mm
- Load capacity: 200 kg; according to IEC 60601 tested with fourfold static load



Design:


2 different variants for artificial leather eligible for upholstery covering:

- Skai® Pandoria Plus with velvet-smooth surface
- Skai® Toronto EN with innovative staynu-technology



Available as an option (see pricelist for additional options)

Mobility

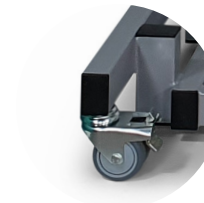

Art. no. 077N
Movable on twin castors
Ø 75 mm, each with total lock



Art. no. 049
Movable on twin castors, Ø 100 mm
with central lock




Art. no. 044N
Wheel-lifting/fixing mechanism
with 4 foot levers for controlled
lowering

Paper roll holder



Art. no. 010F/K
Paper roll holder, mounted under head
part (not protruding), for width 700 mm

Nose opening

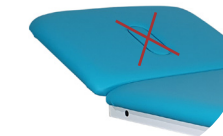

Art. no. 1615-EF
Loss of nose opening and upholstery
part to close

Upholstery options

Art. no. B1770
Soft covering for pleasant resting comfort,
for width 700 mm, thickness 80 mm

Art. no. B1772
Upholstery medium-soft (cold foam cover),
for width 700 mm, thickness 70 mm

