

THERAPY TABLES

2021.1



SERIES 2800XL AND 4800E



INCLINATION
ADJUSTMENT




ROOF POSITION



GYMNASTIC TABLES
(BOBATH & VOJTA)



THERAPY TABLES

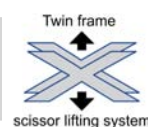



Therapy tables

Content

Good reasons for a therapy table of series 2800XL and 4800E	3
---	---

Therapy tables - Series 2800XL

3-sectioned therapy tables:

• Model 2815XLE	6
• Model 2815XL/H	8
• Model 2820XLE	10
• Model 2820XL/H	12

3-sectioned therapy tables, all upholstery sections adjustable:

• Model 2830XLE	14
• Model 2830XL/H	16
• Model 2890XLE	18
• Model 2890XL/H	20

4-sectioned therapy tables, all upholstery sections adjustable:

• Model 2895XLE	22
• Model 2895XL/H	24
• Model 2850XLE	26
• Model 2850XL/H	28
• Model 2860XLE	30
• Model 2860XL/H	32

6-sectioned therapy tables, all upholstery sections adjustable:

• Model 2880XLE	34
• Model 2880XL/H	36

Therapy tables - Series 4800E

• Model 4815E	38
• Model 4840E	40

Gymnastic tables for therapy acc. to Bobath and Vojta

• Series 2800XLE	42
• Series 4800E	44


Additional equipment	46
----------------------	----

Therapy tables

Good reasons for a therapy table of series 2800XL and 4800E

Electric motorized height adjustment Series 2800XL und 4800E

- Electronical height adjustment by AC-motor with a pressure capacity of 10000 N
- Constant adjustment speed even in case of high load
- EPR-safety power cable with strain relief, length: 4 m
- **Integrated safety devices**
(according to the newest advice of BfArM)
 - **Automatic blocking** of electromotive drive
 - **Active freewheeling downwards**, automatic declutch in case of resistance
- High load capacity
- Series 4800E additionally: High adjustment range

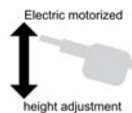

Hydraulic height adjustment Series 2800XL

- Manual foot operated height adjustment by a hydraulic pump with a pressure capacity of 10000 N
- Pedals for foot-controlled height positioning on both sides of the table
- Load capacity: 250 kg



Dewert-hygiene-equipment

Enables optimal cleaning and disinfection



- Hinge lid made of identical upholstery covering
- Artificial leather underneath the upholstery
- Vent holes underneath the upholstery, closed by special caps with valve function
- Easy hygiene because of an open construction: good accessibility of components, no plastic coverings which are breakable in the ageing process
- Option: Innovative artificial leather Skai® Toronto EN with staynu



Therapy tables

Good reasons for a therapy table of series 2800XL and 4800E




Durable ✓



Robust ✓



Easy-care ✓

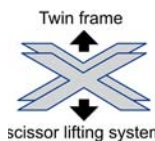

Twin frame scissor lifting system ✓




High load capacity ✓



Universal use ✓



Proven Dewert-quality

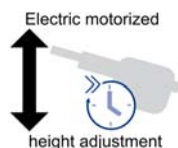

High adjustment speed ✓

Technology:



- Extremely stable twin frame scissor lifting system, with high torsional stiffness due to thick-walled square tubes with big dimensions, also for under and upper frame and scissor

- High stability due to massive 20 mm axes

- Widely-spaced feet and a high deadweight guarantee stability

- Robust frame construction which is proven in demanding everyday life of hospitals for years

- Abrasion-resistant guide system at under and upper frame due to the use of roller bearings



Upholstery:

- Thickness of upholstery: 60 mm
- Multiple glued, high-strength laminated wood panel (thickness: 18 mm) no breaking chipboards in case of pressure on the edges
- Two-layer homogenous polyurethane foam, density 46 and 75, Ökotex-certified
- 2 different variants for artificial leather eligible for upholstery covering:
 - Skai® Pandoria Plus with velvet-smooth surface
 - Skai® Toronto EN with innovative staynu-technology



Cleaning kit for Skai products

(Art. no. 050-1000)

Consists of:

- Spray bottle with 200 ml special Skai cleansing agent
- 1 microfibre care cloth 40 x 40 cm

Therapy tables

Good reasons for a therapy table of series 2800XL and 4800E


... with 3-sectioned lying surface:

- Adjustable head part either positive and negative or only positive adjustable with tubular steel frame for high torsion rigidity
- Pressure stable gas springs ensure working at head part and foot part **without vibrations**
- Operating lever to manipulate the gas spring easily accessible
- **Variable length** due to foldable head part (partly: can be hinged vertical)



... with 3-, 4- or 6-sectioned lying surface, all upholstery sections adjustable:

- Head part positive and negative adjustable
- Adjustment of upholstery sections either by pressure stable gas springs or partly by a self-fixing system of 2 guide rods, without clearance (depending on model)
- Operating levers to manipulate the gas spring easily accessible
- Middle part adjustable to **roof position**, **Trendelenburg position** possible
- **Models 2890 and 2895:**
 - Multifunctional and independent adjustment of upholstery sections:
Roof position, Trendelenburg, reverse Trendelenburg

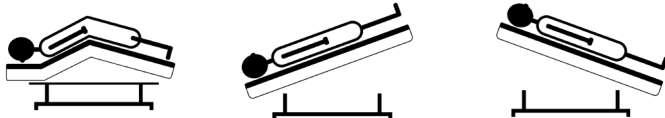

... with 1-sectioned lying surface:

- Gymnastic tables for therapy according to Bobath and Vojta
- Width of lying surface eligible:
100 cm, 110 cm and 120 cm





Series 2800XL

3-sectioned therapy table with electric motorized height adjustment

Model 2815-00XLE white powder coated (RAL 9010)

Model 2815-03XLE white aluminium powder coated (RAL 9006)

Model 2815-04XLE grey aluminium powder coated (RAL 9007)

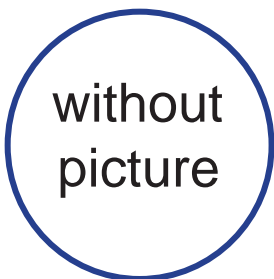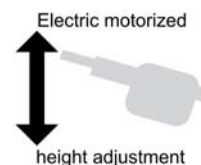

Model 2815-00XL/H (hydraulic height adjustment) with option:
Movable on twin castors with central lock (046)

Basic equipment:

- Lying surface 3-sectioned
- Foot part stepless upwards (+75°) adjustable to sitting position
- Table not movable, including height adjustment screws to compensate for uneven ground
- Electromotive drive control by foot switch (option: hand switch)
- Integrated safety function for automatic blocking of electromotive drive
- Freewheeling downwards, automatic declutch in case of resistance
- Dewert-hygiene-equipment



Under part in case of electric motorized height adjustment

Upholstery:

Lying surface 3-sectioned:

- Head part stepless up- (+45°) and downwards (-35°) adjustable by gas spring
- Length head part: 540 mm, with nose opening and upholstery part to close
- Fixed middle part
- Length middle part: 530 mm
- Foot part stepless upwards (+75°) adjustable to sitting position by gas spring
- Length foot part: 880 mm
- Thickness of upholstery: 60 mm

Technical data:

- Electronical height adjustment: 510 mm to 950 mm in 22 sec. by AC-motor with a pressure capacity of 10000 N
- Width of lying surface: 700 mm
- Length of lying surface: 1950 mm
- Load capacity: 250 kg; according to IEC 60601 tested with fourfold static load



Design:

2 different variants for artificial leather eligible for upholstery covering:

- Skai® Pandoria Plus with velvet-smooth surface
- Skai® Toronto EN with innovative staynu-technology



adjustment

Available as an option (see pricelist for additional options)

Mobility



Art. no. 040N
Movable on comfort twin castors Ø 100 mm, each with total lock, basic height -15 mm



Art. no. 046
Movable on twin castors, Ø 100 mm with central lock, basic height +20 mm



Art. no. B1185
Movable on comfort twin castors, Ø 125 mm with central lock, basic height + 50 mm



Art. no. 045N
Wheel-lifting/fixing mechanism with 4 foot levers for controlled lowering



Additional motor control options



Art. no. 090E
Hand switch (low voltage), with magnet on the back to fasten



Art. no. 091E
Foot switch (low voltage)



Art. no. 092E
1 foot switch fixation mounted at underframe



Art. no. 094E
Foot switch rail, lateral, on both sides, security adjustment




Art. no. 095E
Foot switch rail, surround, security adjustment

Additional safety options



Art. no. B1997E
Foot switch cover, made of brushed VA-steel

Nose opening

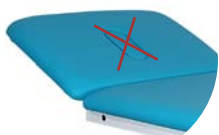

Art. no. 1615-EF
Loss of nose opening and upholstery part to close

Head part options



Art. no. 056
Head part 3-sectioned, arm rests downwards adjustable from -25° to +45° and removable (art. no. 010K not possible)
measures:
Head part: L. 540 mm x W. 370 mm
Arm rests: L. 500 mm x W. 160 mm

Upholstery options

Art. no. 007
Width of upholstery 800 mm

Art. no. B1770
Soft covering for pleasant resting comfort, for width 700 mm, thickness 80 mm

Art. no. B1772
Upholstery medium-soft (cold foam cover), for width 700 mm, thickness 70 mm

Paper roll holder



Art. no. 010K
Paper roll holder, mounted under head part (not protruding), for width 700 mm, not in case of 3-sectioned head part (Art. no. 056)



Art. no. 010K/F
Paper roll holder, mounted at foot part (not protruding), for width 700 mm

Add-on elements



Art. no. B1410
1 piece of side rail, length: 350 mm



Art. no. B1430
1 attachment clamp for side rail, suitable for pipe: Ø 18 mm

Arm rests



Art. no. 088
1 arm rest, fixed, swivelling and removeable, additional revolving and tiltable



Art. no. B1440
1 arm rest for side rail, without attachment clamp, swivelling and removeable, additional revolving and tiltable


Series 2800XL

3-sectioned therapy table with hydraulic height adjustment

Model 2815-00XL/H white powder coated (RAL 9010)

Model 2815-03XL/H white aluminium powder coated (RAL 9006)

Model 2815-04XL/H grey aluminium powder coated (RAL 9007)

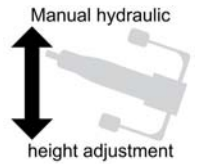

All pictures:
Model 2815-00XL/H
with options:
Movable on twin castors with
central lock (046)



Basic equipment:

- Lying surface 3-sectioned
- Foot part stepless upwards (+75°) adjustable to sitting position
- Table not movable, including height adjustment screws to compensate for uneven ground
- Pedals for foot-controlled height positioning on both sides of the table
- Dewert-hygiene-equipment



Upholstery:

Lying surface 3-sectioned:

- Head part stepless up- (+45°) and downwards (-35°) adjustable by gas spring
- Length head part: 540 mm, with nose opening and upholstery part to close
- Fixed middle part
- Length middle part: 530 mm
- Foot part stepless upwards (+75°) adjustable to sitting position by gas spring
- Length foot part: 880 mm
- Thickness of upholstery: 60 mm



Technical data:


- Manual foot operated height adjustment: 530 mm to 950 mm by a hydraulic pump with a pressure capacity of 10000 N
- Width of lying surface: 700 mm
- Length of lying surface: 1950 mm
- Load capacity: 250 kg; according to IEC 60601 tested with fourfold static load



Design:

2 different variants for artificial leather eligible for upholstery covering:

- Skai® Pandoria Plus with velvet-smooth surface
- Skai® Toronto EN with innovative staynu-technology

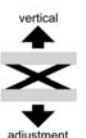

Available as an option (see pricelist for additional options)

Mobility



Art. no. 040N
Movable on comfort twin castors Ø 100 mm, each with total lock, basic height -15 mm



Art. no. 046
Movable on twin castors, Ø 100 mm with central lock, basic height +20 mm



Art. no. B1185
Movable on comfort twin castors, Ø 125 mm with central lock, basic height + 50 mm



Art. no. 045N
Wheel-lifting/fixing mechanism with 4 foot levers for controlled lowering



Paper roll holder



Art. no. 010K
Paper roll holder, mounted under head part (not protruding), for width 700 mm, not in case of 3-sectioned head part (Art. no. 056)



Art. no. 010K/F
Paper roll holder, mounted at foot part (not protruding), for width 700 mm

Upholstery options

Art. no. 007
Width of upholstery 800 mm

Art. no. B1770
Soft covering for pleasant resting comfort, for width 700 mm, thickness 80 mm

Art. no. B1772
Upholstery medium-soft (cold foam cover), for width 700 mm, thickness 70 mm

Head part options



Art. no. 056
Head part 3-sectioned, arm rests downwards adjustable from -25° to +45° and removable (art. no. 010K not possible)
measures:
Head part: L. 540 mm x W. 370 mm
Arm rests: L. 500 mm x W. 160 mm

Nose opening



Art. no. 1615-EF
Loss of nose opening and upholstery part to close

Add-on elements



Art. no. B1410
1 piece of side rail, length: 350 mm



Art. no. B1430
1 attachment clamp for side rail, suitable for pipe: Ø 18 mm

Arm rests

without picture

Art. no. 088
1 arm rest, fixed, swivelling and removeable, additional revolving and tiltable



Art. no. B1440
1 arm rest for side rail, without attachment clamp, swivelling and removeable, additional revolving and tiltable


Series 2800XL

3-sectioned therapy table with electric motorized height adjustment

Model 2820-00XLE white powder coated (RAL 9010)

Model 2820-03XLE white aluminium powder coated (RAL 9006)

Model 2820-04XLE grey aluminium powder coated (RAL 9007)

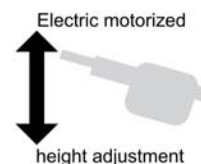

All pictures:
Model 2820-04XLE
with options:
Movable on twin castors with
central lock (046), side rail
(B1410)



Basic equipment:

- Lying surface 3-sectioned
- Head part can be hinged vertical
- Foot part stepless upwards adjustable to sitting position (+80°)
- Table not movable, including height adjustment screws to compensate for uneven ground
- Electromotive drive control by foot switch (option: hand switch)
- Integrated safety function for automatic blocking of electromotive drive
- Freewheeling downwards, automatic declutch in case of resistance
- Dewert-hygiene-equipment



Upholstery:

Lying surface 3-sectioned:

- Head part stepless up- (+30°) and downwards (-90°) adjustable by gas spring
- Length head part: 450 mm, with nose opening and upholstery part to close
- Fixed middle part
- Length middle part: 530 mm
- Foot part stepless upwards (+80°) adjustable to sitting position by gas spring
- Length foot part: 1020 mm
- Thickness of upholstery: 60 mm



Technical data:

- Electronical height adjustment: 510 mm to 950 mm in 22 sec. by AC-motor with a pressure capacity of 10000 N
- Width of lying surface: 700 mm
- Length of lying surface: 2000 mm
- Load capacity: 250 kg; according to IEC 60601 tested with fourfold static load



Design:

2 different variants for artificial leather eligible for upholstery covering:

- Skai® Pandoria Plus with velvet-smooth surface
- Skai® Toronto EN with innovative staynu-technology



Available as an option (see pricelist for additional options)

Mobility



Art. no. 040N
Movable on comfort twin castors Ø 100 mm, each with total lock, basic height -15 mm



Art. no. 046
Movable on twin castors, Ø 100 mm with central lock, basic height +20 mm



Art. no. B1185
Movable on comfort twin castors, Ø 125 mm with central lock, basic height + 50 mm



Art. no. 045N
Wheel-lifting/fixing mechanism with 4 foot levers for controlled lowering



Additional motor control options



Art. no. 090E
Hand switch (low voltage), with magnet on the back to fasten



Art. no. 091E
Foot switch (low voltage)



Art. no. 092E
1 foot switch fixation mounted at underframe

Additional safety options



Art. no. 094E
Foot switch rail, lateral, on both sides, security adjustment



Art. no. B1997E
Foot switch cover, made of brushed VA-steel

Nose opening



Art. no. 1615-EF
Loss of nose opening and upholstery part to close

Paper roll holder



Art. no. 010K
Paper roll holder, mounted under head part (not protruding), for width 700 mm, not in case of 3-sectioned head part (Art. no. 056)



Art. no. 010K/F
Paper roll holder, mounted at foot part (not protruding), for width 700 mm

Upholstery options

Art. no. 007
Width of upholstery 800 mm

Art. no. B1770
Soft covering for pleasant resting comfort, for width 700 mm, thickness 80 mm

Art. no. B1772
Upholstery medium-soft (cold foam cover), for width 700 mm, thickness 70 mm

Add-on elements



Art. no. B1410
1 piece of side rail, length: 350 mm



Art. no. B1430
1 attachment clamp for side rail, suitable for pipe: Ø 18 mm



Art. no. B1440
1 arm rest for side rail, without attachment clamp, swivelling and removeable, additional revolving and tiltable


Series 2800XL

3-sectioned therapy table with hydraulic height adjustment

Model 2820-00XL/H white powder coated (RAL 9010)

Model 2820-03XL/H white aluminium powder coated (RAL 9006)

Model 2820-04XL/H grey aluminium powder coated (RAL 9007)

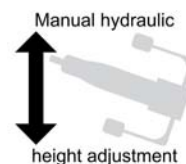

All pictures:
Model 2820-04XL/H
with options:
Movable on twin castors with
central lock (046), paper roll
holder (010K/F)



Basic equipment:

- Lying surface 3-sectioned
- Head part can be hinged vertical
- Foot part stepless upwards adjustable to sitting position (+80°)
- Table not movable, including height adjustment screws to compensate for uneven ground
- Pedals for foot-controlled height positioning on both sides of the table
- Dewert-hygiene-equipment



Upholstery:

Lying surface 3-sectioned:

- Head part stepless up- (+30°) and downwards (-90°) adjustable by gas spring
- Length head part: 450 mm, with nose opening and upholstery part to close
- Fixed middle part
- Length middle part: 530 mm
- Foot part stepless upwards (+80°) adjustable to sitting position by gas spring
- Length foot part: 1020 mm
- Thickness of upholstery: 60 mm



Technical data:

- Manual foot operated height adjustment: 530 mm to 950 mm by a hydraulic pump, pressure capacity of 10000 N
- Width of lying surface: 700 mm
- Length of lying surface: 2000 mm
- Load capacity: 250 kg; according to IEC 60601 tested with fourfold static load



Design:

2 different variants for artificial leather eligible for upholstery covering:

- Skai® Pandoria Plus with velvet-smooth surface
- Skai® Toronto EN with innovative staynu-technology



Available as an option (see pricelist for additional options)

Mobility



Art. no. 040N
Movable on comfort twin castors Ø 100 mm, each with total lock, basic height -15 mm



Art. no. 046
Movable on twin castors, Ø 100 mm with central lock, basic height +20 mm



Art. no. B1185
Movable on comfort twin castors, Ø 125 mm with central lock, basic height + 50 mm



Art. no. 045N
Wheel-lifting/fixing mechanism with 4 foot levers for controlled lowering



Paper roll holder



Art. no. 010K
Paper roll holder, mounted under head part (not protruding), for width 700 mm



Art. no. 010K/F
Paper roll holder, mounted at foot part (not protruding), for width 700 mm

Upholstery options

Art. no. 007
Width of upholstery 800 mm

Art. no. B1770
Soft covering for pleasant resting comfort, for width 700 mm, thickness 80 mm

Art. no. B1772
Upholstery medium-soft (cold foam cover), for width 700 mm, thickness 70 mm

Add-on elements



Art. no. B1410
1 piece of side rail, length: 350 mm



Art. no. B1430
1 attachment clamp for side rail, suitable for pipe: Ø 18 mm

Arm rests



Art. no. B1440
1 arm rest for side rail, without attachment clamp, swivelling and removeable, additional revolving and tiltable

Nose opening



Art. no. 1615-EF
Loss of nose opening and upholstery part to close


Series 2800XL

3-sectioned therapy table with electric motorized height adjustment

Model 2830-00XLE white powder coated (RAL 9010)

Model 2830-03XLE white aluminium powder coated (RAL 9006)

Model 2830-04XLE grey aluminium powder coated (RAL 9007)

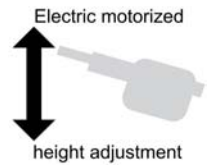

All pictures: Model 2830-00XLE with options: Wheel-lifting/fixing mechanism (045N), surrounding foot switch rail (095E), 3-sectioned head part (056)



Basic equipment:

- Lying surface 3-sectioned
- Middle part stepless adjustable to roof position, Trendelenburg position possible
- Table not movable, including height adjustment screws to compensate for uneven ground
- Electromotive drive control by foot switch (option: hand switch)
- Integrated safety function for automatic blocking of electromotive drive
- Freewheeling downwards, automatic declutch in case of resistance
- Dewert-hygiene-equipment



Upholstery:

Lying surface 3-sectioned:

- Head part stepless up- (+45°) and downwards (-40°) adjustable by gas spring
- Length head part: 540 mm; with nose opening and upholstery part to close
- Middle part upwards manually adjustable to roof position, supported by 2 gas springs, due to 2 guide rods which are fixed by screwing a clamping lever, this adjustment system is without clearance
- Length middle part: 470 mm
- Foot part stepless adjustable (+70°) by a self-fixing system of 2 guide rods, without clearance
- Length foot part: 940 mm
- Thickness of upholstery: 60 mm

Technical data:

- Electronical height adjustment: 510 mm to 950 mm in 22 sec. by AC-motor with a pressure capacity of 10000 N
- Width of lying surface: 700 mm
- Length of lying surface: 1950 mm
- Load capacity: 250 kg; according to IEC 60601 tested with fourfold static load



Design:

2 different variants for artificial leather eligible for upholstery covering:

- Skai® Pandora Plus with velvet-smooth surface
- Skai® Toronto EN with innovative staynu-technology



Available as an option (see pricelist for additional options)

Mobility



Art. no. 040N
Movable on comfort twin castors
Ø 100 mm, each with total lock,
basic height -15 mm



Art. no. 046
Movable on twin castors, Ø 100 mm with
central lock, basic height +20 mm



Art. no. 045N
Wheel-lifting/fixing mechanism
with 4 foot levers for controlled
lowering



Additional motor control options



Art. no. 090E
Hand switch (low voltage),
with magnet on the back to fasten



Art. no. 091E
Foot switch (low voltage)




Art. no. 092E
1 foot switch fixation mounted at
underframe



Art. no. 094E
Foot switch rail, lateral, on both sides,
security adjustment



Art. no. 095E
Foot switch rail, surround,
security adjustment

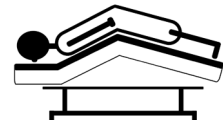

Art. no. 065E
Electric motorized adjustment to roof
position

Additional safety options



Art. no. B1997E
Foot switch cover,
made of brushed VA-steel

Nose opening



Art. no. 1615-EF
Loss of nose opening and upholstery
part to close

Head part options



Art. no. 056
Head part 3-sectioned, arm rests down-
wards adjustable from -25° to +45° and
removable (art. no. 010K not possible)
measures:
Head part: L. 540 mm x W. 370 mm
Arm rests: L. 500 mm x W. 160 mm

Paper roll holder




Art. no. 010K
Paper roll holder, mounted under head
part (not protruding), for width 700 mm,
not in case of 3-sectioned head part
(Art. no. 056)



Art. no. 010K/F
Paper roll holder, mounted at foot part
(not protruding), for width 700 mm

Arm rests

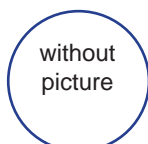

Art. no. 088
1 arm rest, fixed, swivelling and
removeable, additional revolving and
tiltable

Upholstery options

Art. no. 007
Width of upholstery 800 mm

Art. no. B1772
Upholstery medium-soft (cold foam cover),
for width 700 mm, thickness 70 mm


Series 2800XL

3-sectioned therapy table with hydraulic height adjustment

Model 2830-00XL/H white powder coated (RAL 9010)

Model 2830-03XL/H white aluminium powder coated (RAL 9006)

Model 2830-04XL/H grey aluminium powder coated (RAL 9007)

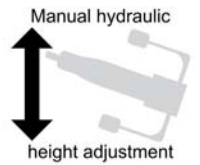

All pictures:
Model 2830-00XL/H
with options:
Wheel-lifting/fixing mechanism
(045N), paper roll holder (010K)



Basic equipment:

- Lying surface 3-sectioned
- Middle part stepless adjustable to roof position, Trendelenburg position possible
- Table not movable, including height adjustment screws to compensate for uneven ground
- Pedals for foot-controlled height positioning on both sides of the table
- Dewert-hygiene-equipment



Upholstery:

Lying surface 3-sectioned:

- Head part stepless up- (+45°) and downwards (-40°) adjustable by gas spring
- Length head part: 540 mm; with nose opening and upholstery part to close
- Middle part upwards manually adjustable to roof position, supported by 2 gas springs, due to 2 guide rods which are fixed by screwing a clamping lever, this adjustment system is without clearance
- Length middle part: 470 mm
- Foot part stepless adjustable (+70°) by a self-fixing system of 2 guide rods, without clearance
- Length foot part: 940 mm
- Thickness of upholstery: 60 mm



Technical data:


- Manual foot operated height adjustment: 530 mm to 950 mm by a hydraulic pump, pressure capacity of 10000 N
- Width of lying surface: 700 mm
- Length of lying surface: 1950 mm
- Load capacity: 250 kg; according to IEC 60601 tested with fourfold static load



Design:

2 different variants for artificial leather eligible for upholstery covering:

- Skai® Pandoria Plus with velvet-smooth surface
- Skai® Toronto EN with innovative staynu-technology

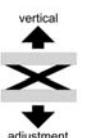

Available as an option (see pricelist for additional options)

Mobility



Art. no. 040N
Movable on comfort twin castors
Ø 100 mm, each with total lock,
basic height -15 mm



Art. no. 046
Movable on twin castors, Ø 100 mm with
central lock, basic height +20 mm



Art. no. 045N
Wheel-lifting/fixing mechanism
with 4 foot levers for controlled
lowering



Paper roll holder



Art. no. 010K
Paper roll holder, mounted under head
part (not protruding), for width 700 mm,
not in case of 3-sectioned head part
(Art. no. 056)



Art. no. 010K/F
Paper roll holder, mounted at foot part
(not protruding), for width 700 mm




Art. no. 056
Head part 3-sectioned, arm rests
downwards adjustable from -25° to
+45° and removable
(art. no. 010K not possible)
measures:
Head part: L. 540 mm x W. 370 mm
Arm rests: L. 500 mm x W. 160 mm

Nose opening



Art. no. 1615-EF
Loss of nose opening and upholstery
part to close

Arm rests

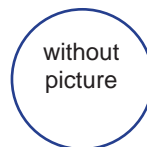

Art. no. 088
1 arm rest, fixed, swivelling and
removeable, additional revolving and
tiltable

Upholstery options

Art. no. 007
Width of upholstery 800 mm

Art. no. B1772
Upholstery medium-soft (cold foam
cover), for width 700 mm,
thickness 70 mm


Series 2800XL

3-sectioned therapy table with electric motorized height adjustment

Model 2890-00XLE white powder coated (RAL 9010)

Model 2890-03XLE white aluminium powder coated (RAL 9006)

Model 2890-04XLE grey aluminium powder coated (RAL 9007)

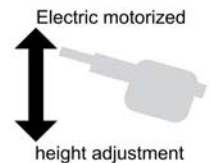

All pictures:
Modell 2890-00XL/H
with option:
Movable on twin castors, each
with lock (077N), paper roll
holder (010F/K);
upholstery colour: 64-713 corn



Basic equipment:

- Lying surface 3-sectioned
- Multifunctional and independent adjustment of upholstery sections
- Possibilities of positioning the surface parts: Roof shape, Trendelenburg, reverse Trendelenburg
- Table not movable, including height adjustment screws to compensate for uneven ground
- Electromotive drive control by foot switch (option: hand switch)
- Integrated safety function for automatic blocking of electromotive drive
- Freewheeling downwards, automatic declutch in case of resistance
- Dewert-hygiene-equipment



Upholstery:

Lying surface 3-sectioned:

- Head part stepless up- (+40°) and downwards (-35°) adjustable by gas spring
- Length head part: 420 mm; with nose opening and upholstery part to close
- Middle part stepless up- (+45°) and downwards (-15°) adjustable by gas spring
- Length middle part: 550 mm
- Foot part stepless up- (+50°) and downwards (-15°) adjustable by gas spring
- Length foot part: 980 mm
- Thickness of upholstery: 60 mm



Technical data:


- Electronical height adjustment: 570 mm to 1010 mm in 22 sec. by AC-motor with a pressure capacity of 10000 N
- Width of lying surface: 700 mm
- Length of lying surface: 1950 mm
- Load capacity: 200 kg; according to IEC 60601 tested with fourfold static load



Design:

2 different variants for artificial leather eligible for upholstery covering:

- Skai® Pandoria Plus with velvet-smooth surface
- Skai® Toronto EN with innovative staynu-technology

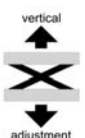

Available as an option (see pricelist for additional options)

Mobility



Art. no. 077N
Movable on twin castors
Ø 75 mm, each with total lock



Art. no. 049
Movable on twin castors, Ø 100 mm
with central lock



Art. no. 044N
Wheel-lifting/fixing mechanism
with 4 foot levers for controlled
lowering



Additional motor control options



Art. no. 090E
Hand switch (low voltage),
with magnet on the back to fasten



Art. no. 091E
Foot switch (low voltage)



Art. no. 092E
1 foot switch fixation mounted at
underframe

Additional safety options



Art. no. B1997E
Foot switch cover,
made of brushed VA-steel



Art. no. 010F/K
Paper roll holder, mounted under head
part (not protruding), for width 700 mm



Art. no. 010F
Paper roll holder, mounted at foot part
(protruding: 70 mm), for width 700 mm

Add-on elements



Art. no. B1410
1 piece of side rail, length: 350 mm



Art. no. B1430
1 attachment clamp for side rail,
suitable for pipe: Ø 18 mm

Arm rests



Art. no. B1440
1 arm rest for side rail, without
attachment clamp, swivelling and
removable, additional revolving and
tiltable

Side guards



Art. no. 1005
1 side guard to fit on side rail, with
attachment clamp, chromium-plated,
length: 800 mm

Nose opening



Art. no. 1615-EF
Loss of nose opening and upholstery
part to close

Upholstery options

Art. no. B1770
Soft covering for pleasant resting comfort, for width 700 mm, thickness 80 mm

Art. no. B1772
Upholstery medium-soft (cold foam cover), for width 700 mm, thickness 70 mm


Series 2800XL

3-sectioned therapy table with hydraulic height adjustment

Model 2890-00XL/H white powder coated (RAL 9010)

Model 2890-03XL/H white aluminium powder coated (RAL 9006)

Model 2890-04XL/H grey aluminium powder coated (RAL 9007)

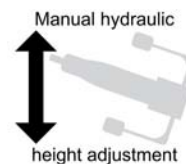

All pictures:
Model 2890-00XL/H
with options:
Movable on twin castors with
central lock (046), paper roll
holder (010K/F)



Basic equipment:

- Lying surface 3-sectioned
- Multifunctional and independent adjustment of upholstery sections
- Possibilities of positioning the surface parts: Roof shape, Trendelenburg, reverse Trendelenburg
- Table not movable, including height adjustment screws to compensate for uneven ground
- Pedals for foot-controlled height positioning on both sides of the table
- Dewert-hygiene-equipment



Upholstery:

Lying surface 3-sectioned:

- Head part stepless up- (+40°) and downwards (-35°) adjustable by gas spring
- Length head part: 420 mm; with nose opening and upholstery part to close
- Middle part stepless up- (+45°) and downwards (-15°) adjustable by gas spring
- Length middle part: 550 mm
- Foot part stepless up- (+50°) and downwards (-15°) adjustable by gas spring
- Length foot part: 980 mm
- Thickness of upholstery: 60 mm

Technical data:

- Manual foot operated height adjustment: 600 mm to 1020 mm by a hydraulic pump, pressure capacity of 10000 N
- Width of lying surface: 700 mm
- Length of lying surface: 1950 mm
- Load capacity: 200 kg; according to IEC 60601 tested with fourfold static load



Design:

2 different variants for artificial leather eligible for upholstery covering:

- Skai® Pandoria Plus with velvet-smooth surface
- Skai® Toronto EN with innovative staynu-technology



Available as an option (see pricelist for additional options)

Mobility



Art. no. 077N
Movable on twin castors
Ø 75 mm, each with total lock



Art. no. 049
Movable on twin castors, Ø 100 mm
with central lock



Art. no. 044N
Wheel-lifting/fixing mechanism
with 4 foot levers for controlled
lowering



Paper roll holder



Art. no. 010F/K
Paper roll holder, mounted under head
part (not protruding), for width 700 mm



Art. no. 010F
Paper roll holder, mounted at foot part
(protruding: 70 mm), for width 700 mm

Nose opening



Art. no. 1615-EF
Loss of nose opening and upholstery
part to close

Add-on elements



Art. no. B1410
1 piece of side rail, length: 350 mm



Art. no. B1430
1 attachment clamp for side rail,
suitable for pipe: Ø 18 mm

Arm rests



Art. no. B1440
1 arm rest for side rail, without
attachment clamp, swivelling and
removable, additional revolving and
tiltable

Side guards



Art. no. 1005
1 side guard to fit on side rail, with
attachment clamp, chromium-plated,
length: 800 mm

Upholstery options

Art. no. B1770
Soft covering for pleasant resting comfort,
for width 700 mm, thickness 80 mm

Art. no. B1772
Upholstery medium-soft (cold foam cover),
for width 700 mm, thickness 70 mm


Series 2800XL

4-sectioned therapy table with electric motorized height adjustment

Model 2895-00XLE white powder coated (RAL 9010)

Model 2895-03XLE white aluminium powder coated (RAL 9006)

Model 2895-04XLE grey aluminium powder coated (RAL 9007)

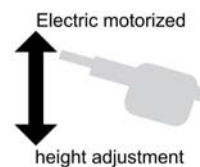

Model 2895-00XLE
with options:
Movable on twin castors with
central lock (046), paper roll
holder (010F/K)



Basic equipment:

- Lying surface 4-sectioned
- Multifunctional and independent adjustment of upholstery sections
- Possibilities of positioning the surface parts: Roof shape, Trendelenburg, reverse Trendelenburg
- Table not movable, including height adjustment screws to compensate for uneven ground
- Electromotive drive control by foot switch (option: hand switch)
- Integrated safety function for automatic blocking of electromotive drive
- Freewheeling downwards, automatic declutch in case of resistance
- Dewert-hygiene-equipment



Upholstery:

Lying surface 4-sectioned:

- Head part stepless up- (+40°) and downwards (-35°) adjustable by gas spring
- Length head part: 420 mm; with nose opening and upholstery part to close
- Middle part stepless up- (+45°) and downwards (-15°) adjustable by gas spring
- Length middle part: 550 mm
- Foot part 2-sectioned, split transversal:
Upper foot segment stepless up- (+50°) and downwards (-15°) adjustable by gas spring
Lower foot segment stepless up- (+35°) and downwards (-30°) adjustable by gas spring
- Length foot part: Upper foot segment: 550 mm | Lower foot segment: 420 mm
- Thickness of upholstery: 60 mm



Technical data:

- Electronical height adjustment: 570 mm to 1010 mm in 22 sec. by AC-motor with a pressure capacity of 10000 N
- Width of lying surface: 700 mm
- Length of lying surface: 1950 mm
- Load capacity: 200 kg; according to IEC 60601 tested with fourfold static load



Design:

2 different variants for artificial leather eligible for upholstery covering:

- Skai® Pandoria Plus with velvet-smooth surface
- Skai® Toronto EN with innovative staynu-technology



Available as an option (see pricelist for additional options)

Mobility



Art. no. 077N
Movable on twin castors
Ø 75 mm, each with total lock



Art. no. 049
Movable on twin castors, Ø 100 mm
with central lock



Art. no. 044N
Wheel-lifting/fixing mechanism
with 4 foot levers for controlled
lowering



Additional motor control options



Art. no. 090E
Hand switch (low voltage),
with magnet on the back to fasten



Art. no. 091E
Foot switch (low voltage)



Art. no. 092E
1 foot switch fixation mounted at
underframe

Additional safety options



Art. no. B1997E
Foot switch cover,
made of brushed VA-steel

Paper roll holder



Art. no. 010F/K
Paper roll holder, mounted under head
part (not protruding), for width 700 mm

Nose opening



Art. no. 1615-EF
Loss of nose opening and upholstery
part to close

Upholstery options

Art. no. B1770
Soft covering for pleasant resting comfort,
for width 700 mm, thickness 80 mm

Art. no. B1772
Upholstery medium-soft (cold foam cover),
for width 700 mm, thickness 70 mm


Series 2800XL

4-sectioned therapy table with hydraulic height adjustment

Model 2895-00XL/H white powder coated (RAL 9010)

Model 2895-03XL/H white aluminium powder coated (RAL 9006)

Model 2895-04XL/H grey aluminium powder coated (RAL 9007)

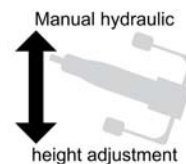

Under part in case of hydraulic height adjustment,
with options:
Movable on twin castors with central lock (046)

without
picture

Model 2895-00XLE (electric motorized height adjustment) with options:
Movable on twin castors with central lock (046),
paper roll holder (010F/K)



Basic equipment:

- Lying surface 4-sectioned
- Multifunctional and independent adjustment of upholstery sections
- Possibilities of positioning the surface parts: Roof shape, Trendelenburg, reverse Trendelenburg
- Table not movable, including height adjustment screws to compensate for uneven ground
- Pedals for foot-controlled height positioning on both sides of the table
- Dewert-hygiene-equipment



Upholstery:

Lying surface 4-sectioned:

- Head part stepless up- (+40°) and downwards (-35°) adjustable by gas spring
- Length head part: 420 mm; with nose opening and upholstery part to close
- Middle part stepless up- (+45°) and downwards (-15°) adjustable by gas spring
- Length middle part: 550 mm
- Foot part 2-sectioned, split transversal:
Upper foot segment stepless up- (+50°) and downwards (-15°) adjustable by gas spring
Lower foot segment stepless up- (+35°) and downwards (-30°) adjustable by gas spring
- Length foot part: Upper foot segment: 550 mm | Lower foot segment: 420 mm
- Thickness of upholstery: 60 mm



Technical data:

- Manual foot operated height adjustment: 600 mm to 1020 mm by a hydraulic pump, pressure capacity of 10000 N
- Width of lying surface: 700 mm
- Length of lying surface: 1950 mm
- Load capacity: 200 kg; according to IEC 60601 tested with fourfold static load



Design:

2 different variants for artificial leather eligible for upholstery covering:

- Skai® Pandoria Plus with velvet-smooth surface
- Skai® Toronto EN with innovative staynu-technology



Available as an option (see pricelist for additional options)

Mobility



Art. no. 077N
Movable on twin castors
Ø 75 mm, each with total lock



Art. no. 049
Movable on twin castors, Ø 100 mm
with central lock



Art. no. 044N
Wheel-lifting/fixing mechanism
with 4 foot levers for controlled
lowering



Paper roll holder



Art. no. 010F/K
Paper roll holder, mounted under head
part (not protruding), for width 700 mm

Nose opening



Art. no. 1615-EF
Loss of nose opening and upholstery
part to close

Upholstery options

Art. no. B1770
Soft covering for pleasant resting comfort,
for width 700 mm, thickness 80 mm

Art. no. B1772
Upholstery medium-soft (cold foam cover),
for width 700 mm, thickness 70 mm



Series 2800XL

4-sectioned therapy table with electric motorized height adjustment

Model 2850-00XLE white powder coated (RAL 9010)

Model 2850-03XLE white aluminium powder coated (RAL 9006)

Model 2850-04XLE grey aluminium powder coated (RAL 9007)

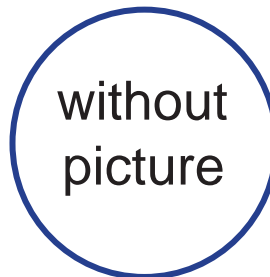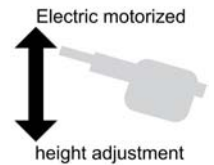

Basic equipment:

- Lying surface 4-sectioned
- Middle part adjustable to roof position, Trendelenburg position possible
- Foot part 2-sectioned, split transversal
- Table not movable, including height adjustment screws to compensate for uneven ground
- Electromotive drive control by foot switch (option: hand switch)
- Integrated safety function for automatic blocking of electromotive drive
- Freewheeling downwards, automatic declutch in case of resistance
- Dewert-hygiene-equipment



Upholstery:

Lying surface 4-sectioned:

- Head part stepless up- (+45°) and downwards (-40°) adjustable by gas spring
- Length head part: 540 mm; with nose opening and upholstery part to close
- Middle part upwards manually adjustable to roof position, supported by 2 gas springs, due to 2 guide rods which are fixed by screwing a clamping lever, this adjustment system is without clearance
- Length middle part: 470 mm
- Foot part 2-sectioned, split transversal
- Complete foot section stepless adjustable (+60°) by a self-fixing system of 2 guide rods, without clearance
- Lower foot section additionally stepless adjustable (-65°) by gas spring (roof position)
- Length foot part: Upper foot section: 405 mm | Lower foot section: 485 mm
- Thickness of upholstery: 60 mm



Technical data:


- Electronical height adjustment: 510 mm to 950 mm in 22 sec. by AC-motor with a pressure capacity of 10000 N
- Width of lying surface: 700 mm
- Length of lying surface: 1950 mm
- Load capacity: 250 kg; according to IEC 60601 tested with fourfold static load



Design:

2 different variants for artificial leather eligible for upholstery covering:

- Skai® Pandoria Plus with velvet-smooth surface
- Skai® Toronto EN with innovative staynu-technology

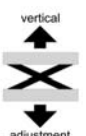

Available as an option (see pricelist for additional options)

Mobility



Art. no. 040N
Movable on comfort twin castors
Ø 100 mm, each with total lock,
basic height -15 mm



Art. no. 046
Movable on twin castors, Ø 100 mm
with central lock, basic height +20 mm



Art. no. 045N
Wheel-lifting/fixing mechanism
with 4 foot levers for controlled
lowering



Additional motor control options



Art. no. 090E
Hand switch (low voltage),
with magnet on the back to fasten



Art. no. 091E
Foot switch (low voltage)



Art. no. 092E
1 foot switch fixation mounted at
underframe



Art. no. 094E
Foot switch rail, lateral, on both sides,
security adjustment



Art. no. 095E
Foot switch rail, surround,
security adjustment

Additional safety options



Art. no. B1997E
Foot switch cover,
made of brushed VA-steel

Paper roll holder



Art. no. 010K
Paper roll holder, mounted under head
part (not protruding), for width 700 mm,
not in case of 3-sectioned head part
(Art. no. 056)

Head part options



Art. no. 056
Head part 3-sectioned, arm rests
downwards adjustable from -25° to
+45° and removable
(art. no. 010K not possible)
measures:
Head part: L. 540 mm x W. 370 mm
Arm rests: L. 500 mm x W. 160 mm

Upholstery options

Art. no. 008
Width of upholstery 800 mm

Art. no. B1772
Upholstery medium-soft (cold foam
cover), for width 700 mm,
thickness 70 mm

Nose opening



Art. no. 1615-EF
Loss of nose opening and upholstery
part to close


Series 2800XL

4-sectioned therapy table with hydraulic height adjustment

Model 2850-00XL/H white powder coated (RAL 9010)

Model 2850-03XL/H white aluminium powder coated (RAL 9006)

Model 2850-04XL/H grey aluminium powder coated (RAL 9007)

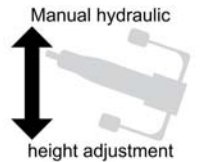

without
picture



Basic equipment:

- Lying surface 4-sectioned
- Middle part adjustable to roof position, Trendelenburg position possible
- Foot part 2-sectioned, split transversal
- Table not movable, including height adjustment screws to compensate for uneven ground
- Pedals for foot-controlled height positioning on both sides of the table
- Dewert-hygiene-equipment



Upholstery:

Lying surface 4-sectioned:

- Head part stepless up- (+45°) and downwards (-40°) adjustable by gas spring
- Length head part: 540 mm; with nose opening and upholstery part to close
- Middle part upwards manually adjustable to roof position, supported by 2 gas springs, due to 2 guide rods which are fixed by screwing a clamping lever, this adjustment system is without clearance
- Length middle part: 470 mm
- Foot part 2-sectioned, split transversal
- Complete foot section stepless adjustable (+60°) by a self-fixing system of 2 guide rods, without clearance
- Lower foot section additionally stepless adjustable (-65°) by gas spring (roof position)
- Length foot part: Upper foot section: 405 mm | Lower foot section: 485 mm
- Thickness of upholstery: 60 mm



Technical data:


- Manual foot operated height adjustment: 530 mm to 950 mm by a hydraulic pump, pressure capacity of 10000 N
- Width of lying surface: 700 mm
- Length of lying surface: 1950 mm
- Load capacity: 250 kg; according to IEC 60601 tested with fourfold static load



Design:

2 different variants for artificial leather eligible for upholstery covering:

- Skai® Pandoria Plus with velvet-smooth surface
- Skai® Toronto EN with innovative staynu-technology

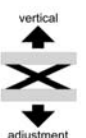

Available as an option (see pricelist for additional options)

Mobility



Art. no. 040N
Movable on comfort twin castors
Ø 100 mm, each with total lock,
basic height -15 mm



Art. no. 046
Movable on twin castors, Ø 100 mm with
central lock, basic height +20 mm



Art. no. 045N
Wheel-lifting/fixing mechanism
with 4 foot levers for controlled
lowering



Paper roll holder



Art. no. 010K
Paper roll holder, mounted under head
part (not protruding), for width 700 mm,
not in case of 3-sectioned head part
(Art. no. 056)

Head part options



Art. no. 056
Head part 3-sectioned, arm rests down-
wards adjustable from -25° to +45° and
removable
(art. no. 010K not possible)
measures:
Head part: L. 540 mm x W. 370 mm
Arm rests: L. 500 mm x W. 160 mm

Upholstery options

Art. no. 008
Width of upholstery 800 mm

Art. no. B1772
Upholstery medium-soft (cold foam
cover), for width 700 mm,
thickness 70 mm

Nose opening



Art. no. 1615-EF
Loss of nose opening and upholstery
part to close


Series 2800XL

4-sectioned therapy table with electric motorized height adjustment

Model 2860-00XLE white powder coated (RAL 9010)

Model 2860-03XLE white aluminium powder coated (RAL 9006)

Model 2860-04XLE grey aluminium powder coated (RAL 9007)

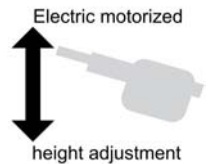

All pictures:
Model 2860-00XLE
in basic equipment



Basic equipment:

- Lying surface 4-sectioned
- Middle part adjustable to roof position, Trendelenburg position possible
- Foot part 2-sectioned, split longitudinal
- Table not movable, including height adjustment screws to compensate for uneven ground
- Electromotive drive control by foot switch (option: hand switch)
- Integrated safety function for automatic blocking of electromotive drive
- Freewheeling downwards, automatic declutch in case of resistance
- Dewert-hygiene-equipment



Upholstery:

Lying surface 4-sectioned:

- Head part stepless up- (+45°) and downwards (-40°) adjustable by gas spring
- Length head part: 540 mm; with nose opening and upholstery part to close
- Middle part upwards manually adjustable to roof position, supported by 2 gas springs, due to 2 guide rods which are fixed by screwing a clamping lever, this adjustment system is without clearance
- Length middle part: 470 mm
- Foot part 2-sectioned, split longitudinal | Length foot part: 940 mm
- Every foot section stepless adjustable (+65°) by a self-fixing system of 2 guide rods, without clearance (foot sections independently adjustable)
- Thickness of upholstery: 60 mm



Technical data:


- Electronical height adjustment: 510 mm to 950 mm in 22 sec. by AC-motor with a pressure capacity of 10000 N
- Width of lying surface: 700 mm
- Length of lying surface: 1950 mm
- Load capacity: 250 kg; according to IEC 60601 tested with fourfold static load



Design:

2 different variants for artificial leather eligible for upholstery covering:

- Skai® Pandoria Plus with velvet-smooth surface
- Skai® Toronto EN with innovative staynu-technology

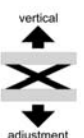

Available as an option (see pricelist for additional options)

Mobility



Art. no. 040N
Movable on comfort twin castors
Ø 100 mm, each with total lock,
basic height -15 mm



Art. no. 046
Movable on twin castors, Ø 100 mm
with central lock, basic height +20 mm



Art. no. 045N
Wheel-lifting/fixing mechanism
with 4 foot levers for controlled
lowering



Additional motor control options



Art. no. 090E
Hand switch (low voltage),
with magnet on the back to fasten



Art. no. 091E
Foot switch (low voltage)



Art. no. 092E
1 foot switch fixation mounted at
underframe

Additional safety options



Art. no. 094E
Foot switch rail, lateral, on both sides,
security adjustment



Art. no. 095E
Foot switch rail, surround,
security adjustment



Art. no. B1997E
Foot switch cover,
made of brushed VA-steel

Paper roll holder



Art. no. 010K
Paper roll holder, mounted under head
part (not protruding), for width 700 mm,
not in case of 3-sectioned head part
(Art. no. 056)

Head part options



Art. no. 056
Head part 3-sectioned, arm rests
downwards adjustable from -25° to
+45° and removable
(art. no. 010K not possible)
measures:
Head part: L. 540 mm x W. 370 mm
Arm rests: L. 500 mm x W. 160 mm

Upholstery options

Art. no. 008
Width of upholstery 800 mm

Art. no. B1772
Upholstery medium-soft (cold foam
cover), for width 700 mm,
thickness 70 mm

Nose opening



Art. no. 1615-EF
Loss of nose opening and upholstery
part to close


Series 2800XL

4-sectioned therapy table with hydraulic height adjustment

Model 2860-00XL/H white powder coated (RAL 9010)

Model 2860-03XL/H white aluminium powder coated (RAL 9006)

Model 2860-04XL/H grey aluminium powder coated (RAL 9007)

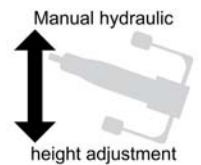

without
picture



Basic equipment:

- Lying surface 4-sectioned
- Middle part adjustable to roof position, Trendelenburg position possible
- Foot part 2-sectioned, split longitudinal
- Table not movable, including height adjustment screws to compensate for uneven ground
- Pedals for foot-controlled height positioning on both sides of the table
- Dewert-hygiene-equipment



Upholstery:

Lying surface 4-sectioned:

- Head part stepless up- (+45°) and downwards (-40°) adjustable by gas spring
- Length head part: 540 mm; with nose opening and upholstery part to close
- Middle part upwards manually adjustable to roof position, supported by 2 gas springs, due to 2 guide rods which are fixed by screwing a clamping lever, this adjustment system is without clearance
- Length middle part: 470 mm
- Foot part 2-sectioned, split longitudinal | Length foot part: 940 mm
- Every foot section stepless adjustable (+65°) by a self-fixing system of 2 guide rods, without clearance (foot sections independently adjustable)
- Thickness of upholstery: 60 mm



Technical data:


- Manual foot operated height adjustment: 530 mm to 950 mm by a hydraulic pump, pressure capacity of 10000 N
- Width of lying surface: 700 mm
- Length of lying surface: 1950 mm
- Load capacity: 250 kg; according to IEC 60601 tested with fourfold static load



Design:

2 different variants for artificial leather eligible for upholstery covering:

- Skai® Pandoria Plus with velvet-smooth surface
- Skai® Toronto EN with innovative staynu-technology

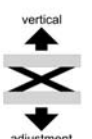

Available as an option (see pricelist for additional options)

Mobility



Art. no. 040N
Movable on comfort twin castors
Ø 100 mm, each with total lock,
basic height -15 mm



Art. no. 046
Movable on twin castors, Ø 100 mm with
central lock, basic height +20 mm



Art. no. 045N
Wheel-lifting/fixing mechanism
with 4 foot levers for controlled
lowering



Paper roll holder



Art. no. 010K
Paper roll holder, mounted under head
part (not protruding), for width 700 mm,
not in case of 3-sectioned head part
(Art. no. 056)

Head part options



Art. no. 056
Head part 3-sectioned, arm rests down-
wards adjustable from -25° to +45° and
removable
(art. no. 010K not possible)
measures:
Head part: L. 540 mm x W. 370 mm
Arm rests: L. 500 mm x W. 160 mm

Upholstery options

Art. no. 008
Width of upholstery 800 mm

Art. no. B1772
Upholstery medium-soft (cold foam
cover), for width 700 mm,
thickness 70 mm

Nose opening



Art. no. 1615-EF
Loss of nose opening and upholstery
part to close


Series 2800XL

6-sectioned therapy table with electric motorized height adjustment

Model 2880-00XLE white powder coated (RAL 9010)

Model 2880-03XLE white aluminium powder coated (RAL 9006)

Model 2880-04XLE grey aluminium powder coated (RAL 9007)

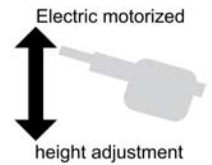

All pictures:
Model 2880-00E
with options:
Wheel-lifting/fixing mechanism (045N), surrounding foot control (095E)



Basic equipment:

- Lying surface 6-sectioned
- Middle part adjustable to roof position, Trendelenburg position possible
- Foot part 4-sectioned, split transversal and longitudinal
- Table not movable, including height adjustment screws to compensate for uneven ground
- Electromotive drive control by foot switch (option: hand switch)
- Integrated safety function for automatic blocking of electromotive drive
- Freewheeling downwards, automatic declutch in case of resistance
- Dewert-hygiene-equipment



Upholstery:

Lying surface 6-sectioned:

- Head part stepless up- (+45°) and downwards (-40°) adjustable by gas spring
- Length head part: 540 mm; with nose opening and upholstery part to close
- Middle part upwards manually adjustable to roof position, supported by 2 gas springs, due to 2 guide rods which are fixed by screwing a clamping lever, this adjustment system is without clearance
- Length middle part: 470 mm
- Foot part 4-sectioned, split transversal and longitudinal
- Every foot section stepless adjustable (+60°) by a self-fixing system of 2 guide rods, without clearance (foot sections independently adjustable)
- Lower foot sections additionally stepless adjustable (-65°) by gas spring
- Length foot part: Upper foot section: 405 mm | Lower foot section: 485 mm
- Thickness of upholstery: 60 mm



Technical data:


- Electronical height adjustment: 510 mm to 950 mm in 22 sec. by AC-motor with a pressure capacity of 10000 N
- Width of lying surface: 700 mm
- Length of lying surface: 1950 mm
- Load capacity: 250 kg; according to IEC 60601 tested with fourfold static load



Design:

2 different variants for artificial leather eligible for upholstery covering:

- Skai® Pandoria Plus with velvet-smooth surface
- Skai® Toronto EN with innovative staynu-technology

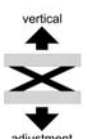

Available as an option (see pricelist for additional options)

Mobility



Art. no. 040N
Movable on comfort twin castors
Ø 100 mm, each with total lock,
basic height -15 mm



Art. no. 046
Movable on twin castors, Ø 100 mm with
central lock, basic height +20 mm



Art. no. 045N
Wheel-lifting/fixing mechanism
with 4 foot levers for controlled
lowering



Additional motor control options



Art. no. 090E
Hand switch (low voltage),
with magnet on the back to fasten



Art. no. 091E
Foot switch (low voltage)




Art. no. 092E
1 foot switch fixation mounted at
underframe



Art. no. 094E
Foot switch rail, lateral, on both sides,
security adjustment



Art. no. 095E
Foot switch rail, surround,
security adjustment

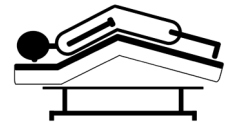

Art. no. 065E
Electric motorized adjustment to roof
position

Additional safety options



Art. no. B1997E
Foot switch cover,
made of brushed VA-steel

Paper roll holder



Art. no. 010K
Paper roll holder, mounted under head
part (not protruding), for width 700 mm,
not in case of 3-sectioned head part
(Art. no. 056)

Nose opening




Art. no. 1615-EF
Loss of nose opening and upholstery
part to close

Head part options



Art. no. 056
Head part 3-sectioned, arm rests down-
wards adjustable from -25° to +45° and
removable
(art. no. 010K not possible)
measures:
Head part: L. 540 mm x W. 370 mm
Arm rests: L. 500 mm x W. 160 mm

Arm rests

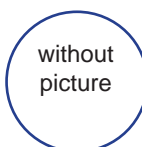

Art. no. 088
1 arm rest, fixed, swivelling and
removeable, additional revolving and
tiltable

Upholstery options

Art. no. 009
Width of upholstery 800 mm

Art. no. B1772
Upholstery medium-soft (cold foam
cover), for width 700 mm,
thickness 70 mm



Series 2800XL

6-sectioned therapy table with hydraulic height adjustment

Model 2880-00XL/H white powder coated (RAL 9010)

Model 2880-03XL/H white aluminium powder coated (RAL 9006)

Model 2880-04XL/H grey aluminium powder coated (RAL 9007)

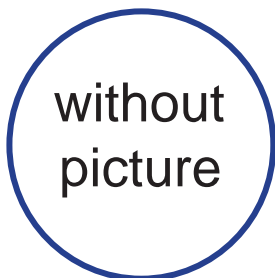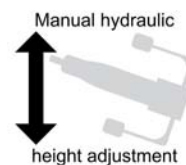

Model 2880-00E (electric motorized height adjustment) with options:
Wheel-lifting/fixing mechanism (045N), surrounding foot control (095E)



Basic equipment:

- Lying surface 6-sectioned
- Middle part adjustable to roof position, Trendelenburg position possible
- Foot part 4-sectioned, split transversal and longitudinal
- Table not movable, including height adjustment screws to compensate for uneven ground
- Pedals for foot-controlled height positioning on both sides of the table
- Dewert-hygiene-equipment



Upholstery:

Lying surface 6-sectioned:

- Head part stepless up- (+45°) and downwards (-40°) adjustable by gas spring
- Length head part: 540 mm; with nose opening and upholstery part to close
- Middle part upwards manually adjustable to roof position, supported by 2 gas springs, due to 2 guide rods which are fixed by screwing a clamping lever, this adjustment system is without clearance
- Length middle part: 470 mm
- Foot part 4-sectioned, split transversal and longitudinal
- Every foot section stepless adjustable (+60°) by a self-fixing system of 2 guide rods, without clearance (foot sections independently adjustable)
- Lower foot sections additionally stepless adjustable (-65°) by gas spring
- Length foot part: Upper foot section: 405 mm | Lower foot section: 485 mm
- Thickness of upholstery: 60 mm



Technical data:

- Manual foot operated height adjustment: 590 mm to 1010 mm by a hydraulic pump, pressure capacity of 10000 N
- Width of lying surface: 700 mm
- Length of lying surface: 1950 mm
- Load capacity: 200 kg; according to IEC 60601 tested with fourfold static load



Design:

2 different variants for artificial leather eligible for upholstery covering:

- Skai® Pandoria Plus with velvet-smooth surface
- Skai® Toronto EN with innovative staynu-technology



Available as an option (see pricelist for additional options)

Mobility



Art. no. 040N
Movable on comfort twin castors
Ø 100 mm, each with total lock,
basic height -15 mm



Art. no. 046
Movable on twin castors, Ø 100 mm
with central lock, basic height +20 mm



Art. no. 045N
Wheel-lifting/fixing mechanism
with 4 foot levers for controlled
lowering



Paper roll holder



Art. no. 010K
Paper roll holder, mounted under head
part (not protruding), for width 700 mm,
not in case of 3-sectioned head part
(Art. no. 056)

Head part options




Art. no. 056
Head part 3-sectioned, arm rests
downwards adjustable from -25° to
+45° and removable
(art. no. 010K not possible)
measures:
Head part: L. 540 mm x W. 370 mm
Arm rests: L. 500 mm x W. 160 mm

Nose opening



Art. no. 1615-EF
Loss of nose opening and upholstery
part to close

Arm rests

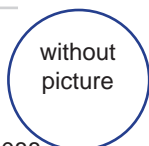

Art. no. B088
1 arm rest, fixed, swivelling and
removeable, additional revolving and
tiltable

Upholstery options

Art. no. 009
Width of upholstery 800 mm

Art. no. B1772
Upholstery medium-soft (cold foam
cover), for width 700 mm,
thickness 70 mm


Series 4800E

3-sectioned therapy table with electric motorized height adjustment

Model 4815-00E white powder coated (RAL 9010)

Model 4815-03E white aluminium powder coated (RAL 9006)

Model 4815-04E grey aluminium powder coated (RAL 9007)

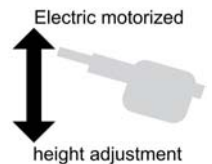

All pictures:
Model 4815-00E
with options:
Movable on twin castors with
central lock (046), surrounding
foot control (095E), paper roll
holder (010F)



Basic equipment:

- Lying surface 3-sectioned
- Foot part stepless adjustable (+75°) to sitting position
- Table not movable, including height adjustment screws to compensate for uneven ground
- High adjustment range
- Electromotive drive control by foot switch (option: hand switch)
- Integrated safety function for automatic blocking of electromotive drive
- Freewheeling downwards, automatic declutch in case of resistance
- Dewert-hygiene-equipment



Table in highest position



Upholstery:

Lying surface 3-sectioned:

- Head part stepless up- (+45°) and downwards (-35°) adjustable by gas spring
- Length head part: 500 mm; with nose opening and upholstery part to close
- Middle part: fixed
- Length middle part: 530 mm
- Foot part stepless upwards adjustable (+75°) to sitting position by gas spring
- Length foot part: 940 mm
- Thickness of upholstery: 60 mm

Table in lowest position



Technical data:


- Electronical height adjustment: 450 mm to 1000 mm in 18 sec. by AC-motor with a pressure capacity of 10000 N
- Width of lying surface: 700 mm
- Length of lying surface: 1970 mm
- Load capacity: 200 kg; according to IEC 60601 tested with fourfold static load



Design:

2 different variants for artificial leather eligible for upholstery covering:

- Skai® Pandora Plus with velvet-smooth surface
- Skai® Toronto EN with innovative staynu-technology

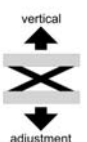

Available as an option (see pricelist for additional options)

Mobility



Art. no. 040N
Movable on comfort twin castors
Ø 100 mm, each with total lock,
basic height -15 mm



Art. no. 046
Movable on twin castors, Ø 100 mm
with central lock, basic height +20 mm



Art. no. 045N
Wheel-lifting/fixing mechanism
with 4 foot levers for controlled
lowering



Additional motor control options



Art. no. 090E
Hand switch (low voltage),
with magnet on the back to fasten



Art. no. 091E
Foot switch (low voltage)



Art. no. 092E
1 foot switch fixation mounted at
underframe

Additional safety options



Art. no. 094E
Foot switch rail, lateral, on both sides,
security adjustment



Art. no. 095E
Foot switch rail, surround,
security adjustment



Art. no. B1997E
Foot switch cover,
made of brushed VA-steel

Paper roll holder



Art. no. 010K
Paper roll holder, mounted under head
part (not protruding), for width 700 mm,
not in case of 3-sectioned head part
(Art. no. 056)




Art. no. 010F
Paper roll holder, mounted at foot part
(protruding: 70 mm), for width 700 mm

Head part options



Art. no. 055
Head part 3-sectioned, arm rests
downwards adjustable from -25° to
+45° and removable
(art. no. 010K not possible)
measures:
Head part: L. 500 mm x W. 370 mm
Arm rests: L. 420 mm x W. 160 mm

Arm rests

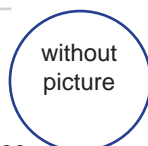

Art. no. 088
1 arm rest, fixed, swivelling and
removeable, additional revolving and
tiltable

Nose opening



Art. no. 1615-EF
Loss of nose opening and upholstery
part to close

Upholstery options

Art. no. 007
Width of upholstery 800 mm

Art. no. B1770
Soft covering for pleasant resting comfort, for width 700 mm, thickness 80 mm

Art. no. B1772
Upholstery medium-soft (cold foam cover), for width 700 mm, thickness 70 mm


Series 4800E

3-sectioned therapy table with electric motorized height adjustment

Model 4840-00E white powder coated (RAL 9010)

Model 4840-03E white aluminium powder coated (RAL 9006)

Model 4840-04E grey aluminium powder coated (RAL 9007)

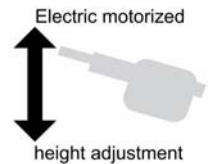

All pictures:
Model 4840-00E
with options:
Wheel-lifting/fixing mechanism
(045N),
surrounding foot control (095E),
head part 3-sectioned (055-10)



Basic equipment:


- Lying surface 3-sectioned
- Multifunctional and independent adjustment of upholstery sections
- Possibilities of positioning the surface parts: Roof shape, Trendelenburg, reverse Trendelenburg
- Table not movable, including height adjustment screws to compensate for uneven ground
- High adjustment range
- Electromotive drive control by foot switch (option: hand switch)
- Integrated safety function for automatic blocking of electromotive drive
- Freewheeling downwards, automatic declutch in case of resistance
- Dewert-hygiene-equipment



Upholstery:

Lying surface 3-sectioned:

- Head part stepless up- (+ 45°) and downwards (-35°) adjustable by gas spring
- Length head part: 450 mm
- Head part with nose opening and upholstery part to close
- Middle part adjustable to roof position by gas spring
- Length middlepart: 530 mm
- Foot part stepless adjustable (+65°) to sitting position by 2 gas springs
- Length foot part: 1020 mm
- Thickness of upholstery: 60 mm

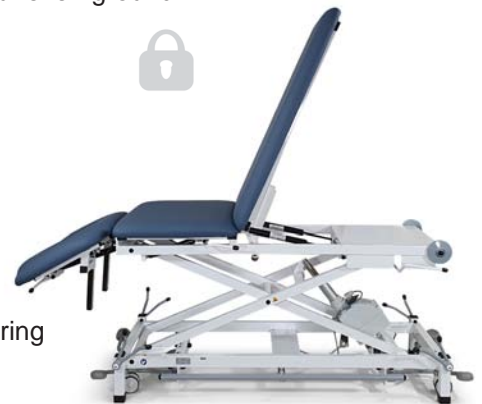

Technical data:


- Electronical height adjustment: 510 mm to 1030 mm in 18 sec. by AC-motor with a pressure capacity of 10000 N
- Width of lying surface: 700 mm
- Length of lying surface: 2000 mm
- Load capacity: 250 kg; according to IEC 60601 tested with fourfold static load



Design:

2 different variants for artificial leather eligible for upholstery covering:

- Skai® Pandoria Plus with velvet-smooth surface
- Skai® Toronto EN with innovative staynu-technology

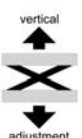

Available as an option (see pricelist for additional options)

Mobility



Art. no. 040N
Movable on comfort twin castors
Ø 100 mm, each with total lock,
basic height -15 mm



Art. no. 046
Movable on twin castors, Ø 100 mm
with central lock, basic height +20 mm



Art. no. 045N
Wheel-lifting/fixing mechanism
with 4 foot levers for controlled
lowering



Additional motor control options



Art. no. 090E
Hand switch (low voltage),
with magnet on the back to fasten



Art. no. 091E
Foot switch (low voltage)



Art. no. 092E
1 foot switch fixation mounted at
underframe

Additional safety options



Art. no. 094E
Foot switch rail, lateral, on both sides,
security adjustment



Art. no. 095E
Foot switch rail, surround,
security adjustment



Art. no. B1997E
Foot switch cover,
made of brushed VA-steel

Paper roll holder



Art. no. 010K
Paper roll holder, mounted under head
part (not protruding), for width 700 mm,
not in case of 3-sectioned head part
(Art. no. 056)




Art. no. 010F
Paper roll holder, mounted at foot part
(protruding: 70 mm), for width 700 mm



Art. no. 055-10
Head part 3-sectioned, arm rests
downwards adjustable from -25° to
+45° and removable
(art. no. 010K not possible)
measures:
Head part: L. 500 mm x W. 370 mm
Arm rests: L. 420 mm x W. 160 mm

Arm rests

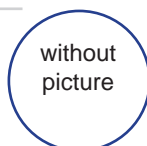

Art. no. B088
1 arm rest, fixed, swivelling and
removeable, additional revolving and
tiltable

Nose opening



Art. no. 1615-EF
Loss of nose opening and upholstery
part to close

Upholstery options

Art. no. 007
Width of upholstery 800 mm

Art. no. B1772
Upholstery medium-soft (cold foam cover), for width 700 mm, thickness 70 mm

Series 2800XL


Gymnastic table for therapy acc. to Bobath and Vojta with electric motorized height adjustment

Model 2800-00XLE* white powder coated (RAL 9010)

Model 2800-03XLE* white aluminium powder coated (RAL 9006)

Model 2800-04XLE* grey aluminium powder coated (RAL 9007)

* Exact name of model depending on width of lying surface (see technical data)

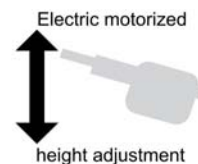

Model 2804-00XLE
with options:
Wheel-lifting/fixing mechanism
(045N), surrounding foot control
(095E), head part stepless
adjustable by gas spring (063)

Basic equipment:

- Lying surface 1-sectioned
- Table not movable, including height adjustment screws to compensate for uneven ground
- Electromotive drive control by foot switch (option: hand switch)
- Integrated safety function for automatic blocking of electromotive drive
- Freewheeling downwards, automatic declutch in case of resistance
- Dewert-hygiene-equipment



Table in lowest position



Both pictures.:
Model 2800-00XLE with
option:
Movable on twin
castors with central lock
(046)

Upholstery:

Lying surface 1-sectioned:

- Length: 1950 mm
- Optional: Head part
- Thickness of upholstery: 60 mm

Technical data:

- Electronical height adjustment: 480 mm to 920 mm in 22 sec. by AC-motor with a pressure capacity of 10000 N
- Width of lying surface: **Model 2800XLE:** 1000 mm | **Model 2802XLE:** 1100 mm | **Model 2804XLE:** 1200 mm
- Length of lying surface: 1950 mm
- Load capacity: 250 kg; according to IEC 60601 tested with fourfold static load



Design:

2 different variants for artificial leather eligible for upholstery covering:

- Skai® Pandoria Plus with velvet-smooth surface
- Skai® Toronto EN with innovative staynu-technology



Available as an option (see pricelist for additional options)

Mobility



Art. no. 040N
Movable on comfort twin castors
Ø 100 mm, each with total lock,
basic height -15 mm



Art. no. 046
Movable on twin castors, Ø 100 mm
with central lock, basic height +20 mm



Art. no. 045N
Wheel-lifting/fixing mechanism
with 4 foot levers for controlled
lowering



Additional motor control options



Art. no. 090E
Hand switch (low voltage),
with magnet on the back to fasten



Art. no. 091E
Foot switch (low voltage)



Art. no. 092E
1 foot switch fixation mounted at
underframe



Art. no. 094E
Foot switch rail, lateral, on both sides,
security adjustment



Art. no. 095E
Foot switch rail, surround,
security adjustment

Additional safety options



Art. no. B1997E
Foot switch cover,
made of brushed VA-steel

Paper roll holder



Art. no. 016K
Paper roll holder, mounted under head
part (not protruding), for width 1000 mm

Head part options



Art. no. 063
Head part, upwards and downwards
adjustable by gas spring

Series 4800E


Gymnastic table for therapy acc. to Bobath and Vojta with electric motorized height adjustment

Model 4800-00E* white powder coated (RAL 9010)

Model 4800-03E* white aluminium powder coated (RAL 9006)

Model 4800-04E* grey aluminium powder coated (RAL 9007)

* Exact name of model depending on width of lying surface (see technical data)

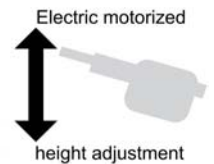

All pictures:
Model 4804-00E
with option:
Wheel-lifting/fixing mechanism
(045N)



Basic equipment:

- Lying surface 1-sectioned
- Table not movable, including height adjustment screws to compensate for uneven ground
- Electromotive drive control by foot switch (option: hand switch)
- High adjustment range
- Integrated safety function for automatic blocking of electromotive drive
- Freewheeling downwards, automatic declutch in case of resistance
- Dewert-hygiene-equipment



Table in lowest position



Table in highest position

Upholstery:

Lying surface 1-sectioned:

- Length: 1950 mm
- Optional: Head part
- Thickness of upholstery: 60 mm

Technical data:

- Electronical height adjustment: 420 mm to 970 mm in 22 sec. by AC-motor with a pressure capacity of 10000 N
- Width of lying surface: **Model 4800E:** 1000 mm | **Model 4802E:** 1100 mm | **Model 4804E:** 1200 mm
- Length of lying surface: 1950 mm
- Load capacity: 250 kg; according to IEC 60601 tested with fourfold static load



Design:

2 different variants for artificial leather eligible for upholstery covering:

- Skai® Pandoria Plus with velvet-smooth surface
- Skai® Toronto EN with innovative staynu-technology



Available as an option (see pricelist for additional options)

Mobility



Art. no. 040N
Movable on comfort twin castors
Ø 100 mm, each with total lock,
basic height -15 mm



Art. no. 046
Movable on twin castors, Ø 100 mm
with central lock, basic height +20 mm



Art. no. 045N
Wheel-lifting/fixing mechanism
with 4 foot levers for controlled
lowering



Additional motor control options



Art. no. 090E
Hand switch (low voltage),
with magnet on the back to fasten



Art. no. 091E
Foot switch (low voltage)



Art. no. 092E
1 foot switch fixation mounted at
underframe



Art. no. 094E
Foot switch rail, lateral, on both sides,
security adjustment



Art. no. 095E
Foot switch rail, surround,
security adjustment

Additional safety options



Art. no. B1997E
Foot switch cover,
made of brushed VA-steel

Paper roll holder



Art. no. 016K
Paper roll holder, mounted under head
part (not protruding), for width 1000 mm

Head part options



Art. no. 063
Head part, upwards and downwards
adjustable by gas spring

K.H. Dewert GmbH
Vollmestraße 7
33649 Bielefeld
Germany

Post office 14 40 28
33634 Bielefeld
Germany

Phone + 49 (0) 5 21 / 4 00 27-0
Fax + 49 (0) 5 21 / 4 00 27-27

Internet: www.khdewert.de
E-Mail: info@khdewert.de
YouTube: K.H. Dewert



05.2021_01

THERAPY TABLES

+ Additional equipment

Cushions



Art. no. 1211-..*
Wedge shaped bolster,
measures (L. x W. x H.):
400 x 400 x 100 - 10 mm

*Add colour number



Art. no. 1113-..*
Half roll bolster,
measures (L. x W. x H.):
570 x 200 x 95 mm



Art. no. 1111-..*
Neck roll,
measures:
L. 570 x Ø 150 mm



Art. no. 1112-..*
Knee roll,
measures:
L. 570 x Ø 200 mm



Art. no. 1314-..*
Cushion, velcro fixation,
ticking with
polyester filling,
measures (L. x W.):
400 x 400 mm

Upholstery protection mat



Art. no. 1890
Protection mat, made of PVC,
thickness: 0.8 mm,
measures: 700 x 600 mm

You will find the whole Dewert product range with additional articles and equipment on our website and in our further catalogues. We look forward to hear from you!

